PRODUCT INFORMATION AND MANUAL

Human GM-CSF Instant ELISA

REF ***BMS283INSTCE***

IVD

Enzyme-linked immunosorbent assay for quantitative
detection of human GM-CSF.

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1 Intended Use

The human GM-CSF Instant ELISA is an enzyme-linked immunosorbent assay for the quantitative detection of human Granulocyte Macrophage Colony Stimulating Factor. **The human GM-CSF Instant ELISA is for in vitro diagnostic use. Not for use in therapeutic procedures.**

2 Summary

Granulocyte - macrophage colony - stimulating factor (GM-CSF) is a small glycoprotein growth factor which stimulates the production and function of neutrophils, eosinophils and monocytes (3) GM-CSF can be produced by a wide variety of tissue types, including fibroblasts, endothelial cells, T-cells, macrophages, mesothelial cells, epithelial cells and many types of tumor cells. In most of these tissues, inflammatory mediators, such as interleukin-1, interleukin-6, tumor necrosis factor or endotoxin, are potent inducers of GM-CSF gene expression, which occurs at least partly by post-transcriptional modification of the GM-CSF mRNA.

The biological effects of GM-CSF are mediated through binding to cell surface receptors, which appear to be widely expressed by hematopoietic cells and also by some non-hematopoietic cells, such as endothelial cells. At least two different functional classes of GM-CSF receptor have been identified. The neutrophil GM-CSF receptor exclusively binds GM-CSF, while interleukin-3 competes for binding of GM-CSF to a second class of receptors detected on some leukemic cell lines (7). Monitoring of GM-CSF has been found relevant in human prostate cancer (14), poorly healing wounds (8), thyroid carcinoma (12), severe mucositis (11), fungal infections (13), AIDS (17), bone marrow transplantation (19, 4, 2, 1), renal cell carcinoma and prostate cancer (18), acute lymphoblastic leukemia (15, 4) pulmonary inflammation (5), hematological malignancies (15, 9), infection (6), lung cancer (5).

3 Principles of the Test

An anti-human GM-CSF coating antibody is adsorbed onto microwells. Human GM-CSF present in the sample or standard binds to antibodies adsorbed to the microwells; a biotin-conjugated anti-human GM-CSF antibody binds to human GM-CSF captured by the first antibody. Streptavidin-HRP binds to the biotin conjugated anti-human GM-CSF.

Following incubation unbound biotin conjugated antihuman GM-CSF and Streptavidin-HRP is removed during a wash step, and substrate solution reactive with HRP is added to the wells.

A coloured product is formed in proportion to the amount of soluble human GM-CSF present in the sample. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from 7 human GM-CSF standard dilutions and human GM-CSF sample concentration determined.

Figure 1

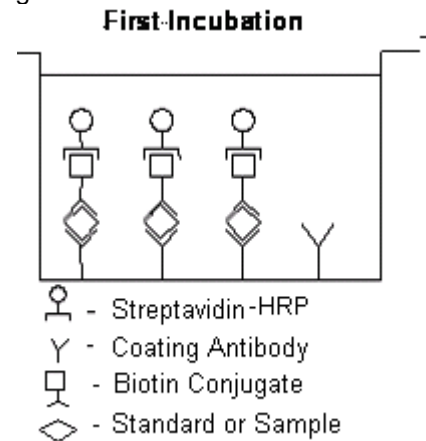
Figure 2


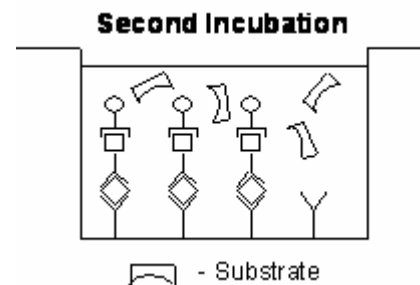
Figure 3

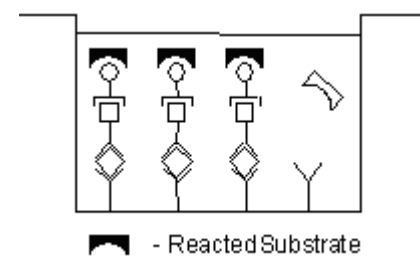4 Reagents Provided

- 1 aluminium pouch with a **Microwell Plate coated** with monoclonal antibody to human GM-CSF, **Biotin-Conjugate** (anti-human GM-CSF monoclonal antibody), **Streptavidin-HRP** and Sample Diluent, lyophilized
- 2 aluminium pouches with a human GM-CSF **Standard curve** (coloured)
- 1 bottle (25 ml) **Wash Buffer Concentrate 20x** (phosphate-buffered saline with 1% Tween 20)
- 1 vial (12 ml) **Sample Diluent** (Use when an external predilution of the samples is needed)
- 1 vial (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 1 vial (15 ml) **Stop Solution** (1M Phosphoric acid)
- 2 **Adhesive Films**

5 Storage Instructions

Store ELISA plate and Standard curves or whole kit at -20°C. The plate and the standard curves can also be removed, stored at -20°C, remaining kit reagents can be stored between 2° and 8°C. Expiry of the kit and reagents is stated on labels.

The expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, the reagent is not contaminated by the first handling.

6 Specimen Collection

Cell culture supernatant, serum, plasma (EDTA) were tested with this assay. Other body fluids might be suitable for use in the assay. Remove serum or plasma from the clot or cells as soon as possible after clotting and separation.

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Pay attention to a possible “**Hook Effect**” due to high sample concentrations (see chapter 11).

Samples must be stored frozen at -20°C to avoid loss of bioactive human GM-CSF. If samples are to be run within 24 hours, they may be stored at 2° to 8°C (for sample stability refer to 13).

Avoid repeated freeze-thaw cycles. Prior to assay, the frozen sample should be brought to room temperature slowly and mixed gently.

7 Materials Required But Not Provided

- 5 ml and 10 ml graduated pipettes
- 5 μ l to 1000 μ l adjustable single channel micropipettes with disposable tips
- adjustable multichannel micropipettes (for volumes between 50 μ l and 500 μ l) with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform linear regression analysis

8 Precautions for Use

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statements(s) for specific advice.
- Reagents are intended for in vitro diagnostic use and are not for use in therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipette by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing substrate reagent.

- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5°C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

9 Preparation of Reagents and Samples

The **buffer concentrate** should be brought to room temperature and diluted before starting the test procedure.

9.1 Wash Buffer

If crystals have formed in the **Wash Buffer Concentrate**, warm it gently until they have completely dissolved.

Pour entire contents (25 ml) of the Wash Buffer Concentrate into a clean 500 ml graduated cylinder. Bring to final volume to 500 ml with glass-distilled or deionized water. Mix gently to avoid foaming. The pH of the final solution should adjust to 7.4.

Transfer to a clean wash bottle and store at 2° to 25°C. Please note that Wash Buffer is stable for 30 days.

10 Test Protocol

- **Use plate immediately after removal from -20 °C!**
 - **Do not wait until pellets have completely dissolved before applying samples - the binding reaction in the standard strips starts immediately after addition of water!**
 - **Do not try to dissolve pellets by pipetting up and down in the wells - some parts of the pellet could stick to the tip creating high variation of results.**
 - **Perform the washing step with at least 400 µl of washing buffer as stated in the manual or fill the wells completely - otherwise any pellet residues sticking to the rim of the well will not be removed and create high variation of results.**
 - **Allow the washing buffer to sit in the wells for a few seconds before aspiration.**
 - **Remove covers of the standard strips carefully so that all the lyophilized pellets remain in the wells.**
- a. Determine the number of Microwell Strips required to test the desired number of samples plus Microwell Strips for blanks and standards (coloured). Each sample, standard and blank should be assayed in duplicate. Remove extra Microwell Strips from holder and store in foil bag with the desiccant provided at -20 °C sealed tightly. Place Microwell Strips containing the standard curve in position A1/A2 to H1/H2 (see Table 1).
 - b. Add 100 µl of **distilled water** to the **sample wells**.
 - c. Add **distilled water** to all **standard and blank wells** as indicated on the label of the standard strips (A1, A2 to H1, H2).

Table 1

Table depicting an example of the arrangement of blanks, standards and samples in the microwell strips:

	1	2	3	4
A	Standard 1 (500 pg/ml)	Standard 1 (500 pg/ml)	Sample 1	Sample 1
B	Standard 2 (250 pg/ml)	Standard 2 (250 pg/ml)	Sample 2	Sample 2
C	Standard 3 (125 pg/ml)	Standard 3 (125 pg/ml)	Sample 3	Sample 3
D	Standard 4 (62.5 pg/ml)	Standard 4 (62.5 pg/ml)	Sample 4	Sample 4
E	Standard 5 (31.3 pg/ml)	Standard 5 (31.3 pg/ml)	Sample 5	Sample 5
F	Standard 6 (15.6 pg/ml)	Standard 6 (15.6 pg/ml)	Sample 6	Sample 6
G	Standard 7 (7.8 pg/ml)	Standard 7 (7.8 pg/ml)	Sample 7	Sample 7
H	Blank	Blank	Sample 8	Sample 8

- d. Add 50 µl of each **sample**, in duplicate, to the **designated wells** and mix the contents.
- e. Cover with an adhesive film and incubate at room temperature (18°C to 25°C) for 3 hours, if available on a microplate shaker at 200 rpm.
- f. Remove adhesive film and empty wells. **Wash** the microwell strips **6** times with approximately **400 µl** Wash Buffer per well with thorough aspiration of microwell contents between washes. Allow the Wash Buffer to sit in the wells for about **10 – 15 seconds** before aspiration. Take care not to scratch the surface of the microwells.

After the last wash, tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing or place upside down on a wet absorbent paper for no longer than 15 minutes. Do not allow wells to dry.

- g. Pipette 100 µl of **TMB Substrate Solution** to all wells, including the blank wells.
- h. Incubate the microwell strips at room temperature (18° to 25°C) for about 10 min. Avoid direct exposure to intense light.

The colour development on the plate should be monitored and the substrate reaction stopped (see next point of this protocol) before positive wells are no longer properly recordable. Determination of the ideal time period for colour development has to be done individually for each assay.

It is recommended to add the Stop Solution when the highest standard has developed a dark blue colour. Alternatively the colour development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.6 – 0.65.

- i. Stop the enzyme reaction by quickly pipetting 100 µl of **Stop Solution** into each well, including the blank wells. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2 - 8°C in the dark.

- j. Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the human GM-CSF standards.

Note: In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.

11 Calculation of Results

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20 per cent of the mean.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the human GM-CSF concentration on the abscissa. Draw a best fit curve through the points of the graph (a 5-parameter curve fit is recommended).
- To determine the concentration of circulating human GM-CSF for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding human GM-CSF concentration.
- ***Samples have been diluted 1:2, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 2).**
- **Calculation of samples with a concentration exceeding standard 1 may result in incorrect, low human GM-CSF levels (Hook Effect). Such samples require further external predilution according to expected human GM-CSF values with Sample Diluent in order to precisely quantitate the actual human GM-CSF level.**
- It is suggested that each testing facility establishes a control sample of known human GM-CSF concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.
- A representative standard curve is shown in Figure 4. This curve cannot be used to derive test results. Every laboratory must prepare a standard curve for each group of microwell strips assayed.

* N.B: There is a common dilution factor for samples due to the conjugate which must then be included in the calculation. The samples contribute 100 µl to the final volume per well. These 100 µl are composed of 50 µl of Sample Diluent plus 50 µl of the sample. This is a 1:2 dilution.

The remaining 50 µl to give 150 µl are due to the addition of 50 µl conjugate to all wells. 50 µl Sample Diluent and 50 µl conjugate results in 100 µl reconstitution volume, addition of 50 µl sample (50 µl + 50 µl = 1:2 dilution)

Figure 4

Representative standard curve for human GM-CSF Instant ELISA. Human GM-CSF was diluted in serial 2-fold steps in Sample Diluent. Each symbol represents the mean of 3 parallel titrations. Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

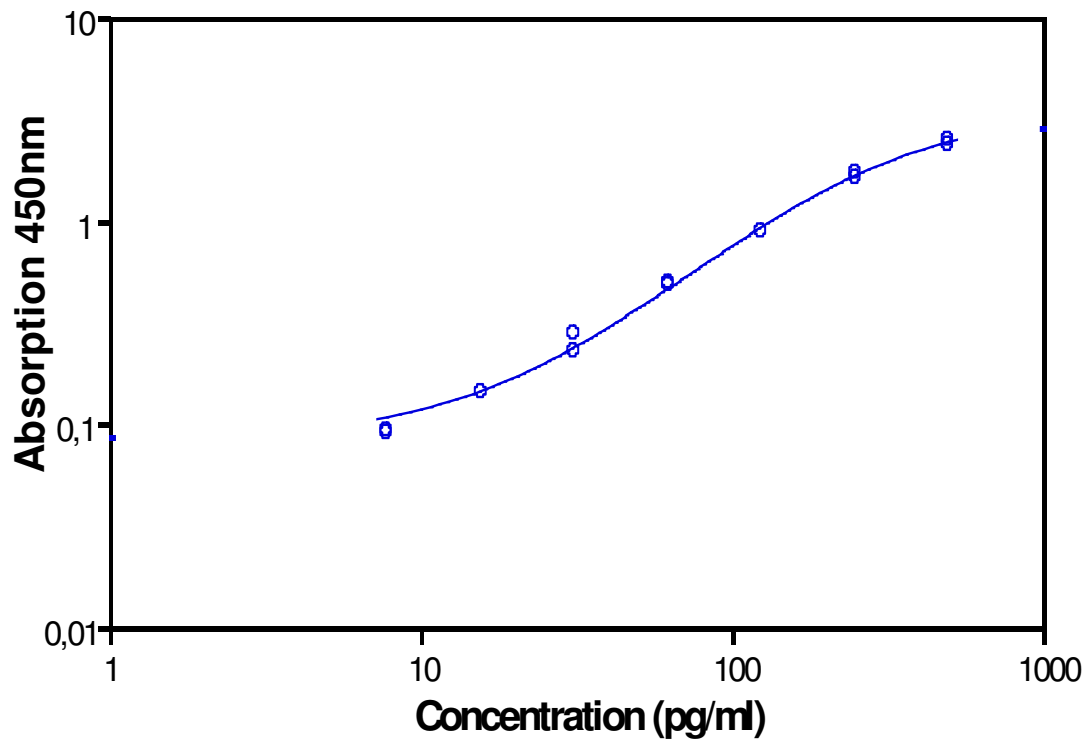
Table 2

Typical data using the human GM-CSF INSTANT ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	human GM-CSF Concentration (pg/ml)	O.D. (450 nm)	O.D. Mean	C.V. (%)
1	500	2.468	2.95	1.8
	500	2.383		
2	250	1.641	1.652	2.5
	250	1.725		
3	125	0.883	0.847	0.6
	125	0.872		
4	62.5	0.501	0.459	2.3
	62.5	0.478		
5	31.3	0.226	0.220	9.8
	31.3	0.275		
6	15.6	0.142	0.111	0.4
	15.6	0.141		
7	7.8	0.092	0.061	0.5
	7.8	0.091		
Blank	0	0.033	0.031	6.5
	0	0.029		

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus colour intensity. Values measured are still valid.

12 Limitations

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.
- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.
- The use of radioimmunotherapy has significantly increased the number of patients with human anti-mouse IgG antibodies (HAMA). HAMA may interfere with assays utilizing murine monoclonal antibodies leading to both false positive and false negative results. Serum samples containing antibodies to murine immunoglobulins can still be analysed in such assays when murine immunoglobulins (serum, ascitic fluid, or monoclonal antibodies of irrelevant specificity) are added to the sample.

13 Performance Characteristics

13.1 Sensitivity

The limit of detection of human GM-CSF defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be 0.6 pg/ml (mean of 6 independent assays).

13.2 Reproducibility

13.2.1 Intra-assay

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 serum samples containing different concentrations of human GM-CSF. 2 standard curves were run on each plate. Data below show the mean human GM-CSF concentration and the coefficient of variation for each sample (see Table 3). The calculated overall intra-assay coefficient of variation was 5,9%.

Table 3

The mean human GM-CSF concentration and the coefficient of variation for each sample.

Positive Sample	Experiment	human GM-CSF Concentration (pg/ml)	Coefficient of Variation (%)
1	1	468	4
	2	368	6
	3	359	9
2	1	371	7
	2	310	6
	3	301	4
3	1	362	6
	2	286	8
	3	272	9
4	1	202	3
	2	184	3
	3	165	9
5	1	146	5
	2	113	6
	3	115	7
6	1	252	5
	2	228	4
	3	192	6
7	1	178	4
	2	140	3
	3	137	9
8	1	144	6
	2	133	4
	3	113	9

13.2.2 Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 serum samples containing different concentrations of human GM-CSF. 2 standard curves were run on each plate. Data below show the mean human GM-CSF concentration and the coefficient of variation calculated on 18 determinations of each sample (see Table 4). The calculated overall inter-assay coefficient of variation was 13.5%.

Table 4

The mean human GM-CSF concentration and the coefficient of variation of each sample

Sample	Mean human GM-CSF Concentration (pg/ml)	Coefficient of Variation (%)
1	399	15.1
2	327	11.6
3	307	15.7
4	184	10.1
5	124	15.0
6	224	13.7
7	152	14.9
8	130	12.2

13.3 Spike Recovery

The spike recovery was evaluated by spiking 4 levels of human GM-CSF into serum. Recoveries were determined in 3 independent experiments with 6 replicates each. The unspiked serum was used as blank in these experiments. Average recovery ranged from 81% to 94% with an overall mean recovery of 88%.

13.4 Dilution Parallelism

Serum samples with different levels of human GM-CSF were analysed at serial 2 fold dilutions with 4 replicates each. The recovery ranged between 79% and 118% with an overall recovery of 102% (see Table 5).

Table 5

Sample	Dilution	Mean human GM-CSF Concentration (pg/ml)		% Recovery of Exp. Val
		Expected Value	Observed Value	
1	1:2	--	846	--
	1:4	423	476	113
	1:8	238	282	118
	1:16	141	112	79
2	1:2	--	796	--
	1:4	398	440	111
	1:8	220	260	118
	1:16	130	112	86
3	1:2	--	522	--
	1:4	261	285	109
	1:8	143	148	104
	1:16	74	63	85
4	1:2	--	492	--
	1:4	246	258	105
	1:8	129	139	108
	1:16	69	61	87

13.5 Sample Stability

13.5.1 Freeze-Thaw Stability

Aliquots of serum samples (unspiked or spiked) were stored at -20 °C and thawed 5 times, and the human GM-CSF levels determined. There was no significant loss of human GM-CSF immunoreactivity detected by freezing and thawing.

13.5.2 Storage Stability

Aliquots of serum samples (spiked or unspiked) were stored at -20 °C, 2-8 °C, room temperature (RT) and at 37 °C, and the human GM-CSF level determined after 24 h.

There was no significant loss of human GM-CSF immunoreactivity detected during storage under above conditions.

13.6 Specificity

To define the specificity of this ELISA several proteins were tested for cross reactivity. There was no cross reactivity observed.

13.7 Expected Values

A panel of sera samples from randomly selected apparently healthy donors (males and females) was tested for human GM-CSF.

There were no detectable human GM-CSF levels found.

Elevated human GM-CSF levels depend on the type of immunological disorder.

14 Bibliography

- 1) Antman KH. GM-CSF in marrow transplantation. *Eur J Cancer*. 1993;29A Suppl 3:S2-6. Review.
- 2) Bolwell BJ. The use of G-CSF and GM-CSF in bone marrow transplantation. *Cleve Clin J Med*. 1993 Jul-Aug;60(4):291-302. Review.
- 3) Burgess AW, Metcalf D. The nature and action of granulocyte-macrophage colony stimulating factors. *Blood* 1980 Dec;56(6):947-58
- 4) Croockewit S. GM-CSF in haematopoietic stem cell transplantation. *Eur J Cancer*. 1999 Aug;35 Suppl 3:S11-3. Review. No abstract available.
- 5) Drings P, Fischer JR. Biology and clinical use of GM-CSF in lung cancer. *Lung*. 1990;168 Suppl:1059-68. Review.
- 6) Freund M, Kleine HD. The role of GM-CSF in infection. *Infection*. 1992;20 Suppl 2:S84-92 Review.
- 7) Griffin JD, Cannistra SA, Sullivan R, Demetri GD, Ernst TJ, Kanakura Y. The biology of GM-CSF: regulation of production and interaction with its receptor. *Int J Cell Cloning* 1990 Jan;8 Suppl 1:35-44; discussion 44-5.
- 8) Groves RW, Schmidt-Lucke JA. Recombinant human GM-CSF in the treatment of poorly healing wounds. *Adv Skin Wound Care*. 2000 May-Jun;13(3 Pt 1):107-12. Review.
- 9) Harmenberg J, Hoglund M, Hellstrom-Lindberg E. G- and GM-CSF in oncology and oncological haematology. *Eur J Haematol Suppl*. 1994;55:1-28. Review.
- 10) Harousseau JL, Wu D. The use of GM-CSF and G-CSF in the treatment of acute leukemias. *Leuk Lymphoma*. 1995 Aug;18(5-6):405-12. Review.
- 11) Hejna M, Brodowicz T, Zielinski CC. Local use of GM-CSF for severe mucositis. *Eur J Cancer*. 1999 Aug;35 Suppl 3:S14-7. Review. No abstract available.

- 12) Hoshi S, Yoshizawa A, Arioka H, Kobayashi N, Kudo K, Niino H. Anaplastic thyroid carcinoma with lung metastasis producing CA 19-9 and GM-CSF. *Nihon Kokyuki Gakkai Zasshi*. 2000 May;38(5):391-7. Review. Japanese.
- 13) Jones TC. Use of granulocyte-macrophage colony stimulating factor (GM-CSF) in prevention and treatment of fungal infections. *Eur J Cancer*. 1999 Aug;35 Suppl 3:S8-10. Review. No abstract available.
- 14) Lim M, Simons JW. Emerging concepts in GM-CSF gene-transduced tumor vaccines for human prostate cancer. *Curr Opin Mol Ther*. 1999 Feb;1(1):64-71. Review. No abstract available.
- 15) Lowenberg B, de Greef GE, Wielenga JJ. Granulocyte-macrophage colony stimulating factor (GM-CSF) in the treatment of hermatological malignancies. *Anticancer Drugs*. 1993 May;4 Suppl 1:17-20. Review.
- 16) Ottmann OG, Hoelzer D. Do G-CSF and GM-CSF contribute to the management of acute lymphoblastic leukemia? *Leukemia*. 1996 Jul;10(7):1111-6. Review. No abstract available.
- 17) Ross SD, DiGeorge A, Connelly JE, Whiting GW, McDonnell N. Safety of GM-CSF in patients with AIDS: a review of the literature. *Pharmacotherapy*. 1998 Nov-Dec;18(6):1290-7 Review.
- 18) Simons JW. Bioactivity of human GM-CSF gene therapy in metastatic renal cell carcinoma and prostate cancer. *Hinyokika Kyo*. 1997 Nov;43(11):821-2. Review. No abstract available.
- 19) Toren A, Nagler A. The implications of granulocyte-monocyte colony-stimulating factor (GM-CSF) in cytotoxicity of bone marrow transplantation. *Cytokines Cell Mol Ther*. 1998 Sep;4(3):199-206. Review.
- 20) Xing Z, Braciak T, Ohkawara Y, Sallenave JM, Foley R, Sime PJ, Jordana M, Graham FL, Gauldie J. Gene transfer for cytokine functional studies in the lung: the multifunctional role of GM-CSF in pulmonary inflammation. *J Leukoc Biol*. 1996 Apr;59(4):481-8. Review.

For literature update refer to **www.bendermedsystems.com**

15 Ordering Information

For orders please contact:

Europe-Headquarters

Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe
phone: +43 1 796 40 40 ext. 114
fax: +43 1 796 40 40 ext. 400
e-mail: order@bendermedsystems.com

USA

Bender MedSystems, Inc.
849 Hinckley Road
Burlingame, CA 94010, USA
toll-free phone: +1 (866) 952 2112
toll-free fax: +1 (877) 952 2112
phone: +1 (650) 952 2112
fax: +1 (650) 952 2252

For technical information please contact:

e-mail: techserv@bendermedsystems.com
www.bendermedsystems.com

Cat.No. BMS283INSTCE human GM-CSF INSTANT ELISA

16 Reagent Preparation Summary

16.1 Wash Buffer

Add **Wash Buffer Concentrate** 20 x (25 ml) to 475 ml distilled water

17 Test Protocol Summary

- Place standard strips in position A1/A2 to H1/H2.
- Add 100 µl distilled water to sample wells.
- Add distilled water, in duplicate, to all standard and blank wells as indicated on the label of the standard strips.
- Add 50 µl sample to designated wells.
- Cover microwell strips and incubate 3 hours at room temperature (18° to 25 °C) if available on a microplate shaker at 200 rpm.
- Empty and wash microwell strips **6** times with 400 µl Wash Buffer.
- Add 100 µl of TMB Substrate Solution to all wells including blank wells.
- Incubate the microwell strips for about 10 minutes at room temperature (18° to 25 °C).
- Add 100 µl Stop Solution to all wells including blank wells.
- Blank microwell reader and measure colour intensity at 450 nm.

Note: Samples have been diluted 1:2, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 2).

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