

PRODUCT INFORMATION AND MANUAL

human sCD27 Instant ELISA

BMS286INST

Enzyme-linked immunosorbent assay for quantitative
detection of human sCD27.

For research use only.

Not for diagnostic or therapeutic procedures.

128 Tests

**human sCD27
BMS286INST**



Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe

TABLE OF CONTENTS

1	Intended Use	3
2	Summary	3
3	Principles of the Test	4
4	Reagents Provided	5
5	Storage Instructions	5
6	Specimen Collection	6
7	Materials Required But Not Provided	7
8	Precautions for Use	8
9	Preparation of Reagents and Samples	10
10	Test Protocol	11
11	Calculation of Results	15
12	Limitations	18
13	Performance Characteristics	19
14	Bibliography	24
15	Ordering Information	26
16	Reagent Preparation Summary	27
17	Test Protocol Summary	28

1 Intended Use

The human sCD27 Instant ELISA is an enzyme-linked immunosorbent assay for the quantitative detection of human sCD27 in cell culture supernatants, human serum, plasma or other body fluids. The human sCD27 Instant ELISA is **for research use only. Not for diagnostic or therapeutic procedures.**

2 Summary

CD27 is a member of the TNF-receptor family (6). The expression of CD27 is restricted to cells of the lymphoid lineage and appears to be related to functional differentiation programs of T- and B-cells (4, 5, 7, 12). Constitutively expressed on T-lymphocytes, CD27 is a costimulatory molecule for a regulatory subset. Induced on B lymphocytes after antigenic challenge, it is a marker of memory cells (9). Following antigen activation in germinal centers, B-cells develop into memory B-cells or plasma cells. CD27 interacts with its ligand CD70 to control the entry of B cells into their pathway (1, 2, 13). Structurally CD27 being a member of the TNF receptor superfamily is related to Fas, CD40 and RANK. These receptors are characterized by an extra cellular domain comprised of two to six repeats of cysteine rich motifs and the existence of a number of structurally related decoy receptors (8).

A soluble isoform of CD27, sCD27 has been described which is found in body fluids and can be used to monitor local and systemic immune activation. In addition elevated serum concentration of sCD27 are found in patients with B-cell malignancies and levels of sCD27 strongly correlate with tumor load (11) and neoplastic B-cell growth (10). It was furthermore shown to be a marker increased in patients suffering from a variety of other immunological disorders such as chronic viral infections like HIV and HTLV-I infection (14), Sjogren's syndrome (3) and thyroid disorders (15).

3 Principles of the Test

An anti-human sCD27 monoclonal coating antibody is adsorbed onto microwells. Human sCD27 present in the sample or standard binds to antibodies adsorbed to the microwells; a biotin-conjugated monoclonal anti-human sCD27 antibody binds to human sCD27 captured by the first antibody. Streptavidin-HRP binds to the biotin conjugated anti-human sCD27.

Following incubation unbound biotin conjugated anti human sCD27 and Streptavidin-HRP is removed during a wash step, and substrate solution reactive with HRP is added to the wells.

A coloured product is formed in proportion to the amount of soluble human sCD27 present in the sample. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from seven human sCD27 standard dilutions and human sCD27 sample concentration determined.

Figure 1

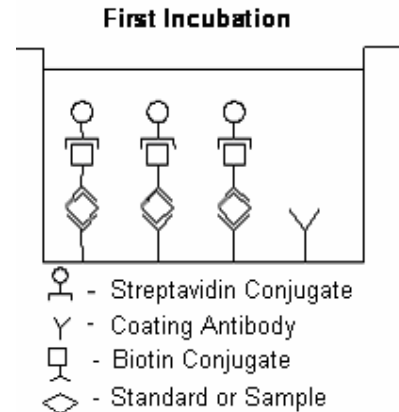


Figure 2


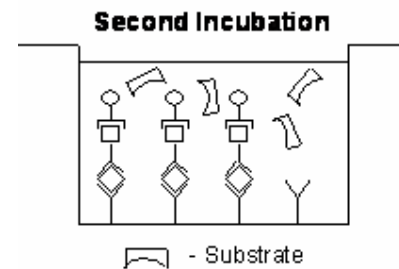


Figure 3

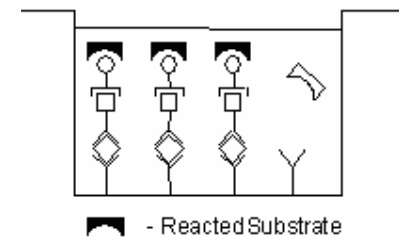

4 Reagents Provided

- 1 aluminium pouch with a **Microwell Plate coated with Monoclonal Antibody** (murine) to human sCD27, **Biotin-Conjugate** (anti-sCD27 monoclonal antibody), Sample Diluent and **Streptavidin-HRP**, lyophilized
- 2 aluminium pouches with a **human sCD27 Standard curve (coloured)**
- 1 bottle (25 ml) **Wash Buffer Concentrate** 20x (phosphate-buffered saline with 1% Tween 20)
- 1 vial (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 1 bottle (50 ml) **Sample Diluent** (Use when an external predilution of the samples is needed)
- 1 vial (15 ml) **Stop Solution** (1M Phosphoric acid)
- 2 adhesive **Plate Covers**

5 Storage Instructions

Store ELISA plate and Standard curves or whole kit at -20°C. The plate and the standard curves can also be removed, stored at -20°C, remaining kit reagents can be stored between 2° and 8°C. Expiry of the kit and reagents is stated on labels.

The expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, the reagent is not contaminated by the first handling.

6 Specimen Collection

Cell culture supernatants, human serum, plasma or other biological samples are suitable for use in the assay. Remove the serum or plasma from the clot or red cells as soon as possible after clotting and separation.

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples must be stored frozen at -20°C to avoid loss of bioactive human sCD27. If samples are to be run within 24 hours, they may be stored at 2° to 8°C (for sample stability refer to 13).

Avoid repeated freeze-thaw cycles. Prior to assay, frozen serum or plasma should be brought to room temperature slowly and mixed gently.

7 Materials Required But Not Provided

- 5 ml and 10 ml graduated pipettes
- 5 μ l to 1000 μ l adjustable single channel micropipettes with disposable tips
- 50 μ l to 300 μ l adjustable multichannel micropipette with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform linear regression analysis

8 Precautions for Use

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statements(s) for specific advice.
- Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipette by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing substrate reagent.

- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5 °C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

9 Preparation of Reagents and Samples

9.1 Wash Buffer

If crystals have formed in the Wash Buffer Concentrate, warm it gently until they have completely dissolved.

Pour entire contents (25 ml) of the Wash Buffer Concentrate into a clean 500 ml graduated cylinder. Bring to final volume to 500 ml with glass-distilled or deionized water. Mix gently to avoid foaming. The pH of the final solution should adjust to 7.4.

Transfer to a clean wash bottle and store at 2° to 25°C. Please note that Wash Buffer is stable for 30 days.

10 Test Protocol

- **Use plate immediately after removal from -20°C!**
- **Do not wait until pellets have completely dissolved before applying samples - the binding reaction in the standard strips starts immediately after addition of water!**
- **Do not try to dissolve pellets by pipetting up and down in the wells - some parts of the pellet could stick to the tip creating high variation of results**
- **Perform the washing step with at least 400 µl of washing buffer as stated in the manual or fill the wells completely - otherwise any pellet residues sticking to the rim of the well will not be removed and create high variation of results**
- **Allow the washing buffer to sit in the wells for a few seconds before aspiration**
- **Remove covers of the standard strips carefully in order that all the lyophilised pellets remain in the wells**
-
- a. Prepare your samples before starting with the test procedure. Dilute serum or plasma samples 1:25 with Sample Diluent according to the following dilution scheme:
10 µl Sample + 240 µl Sample Diluent
- b. Determine the number of microwell Strips required to test the desired number of samples plus microwell Strips for blanks and standards (coloured). Each sample, standard, blank, and optional control sample should be assayed in duplicate. Remove extra microwell Strips from holder and store in foil bag with the desiccant provided at -20°C sealed tightly. Place microwell strips containing the standard curve in position A1/A2 to H1/H2 (see Table 1).
- c. Add 100 µl of distilled water to the sample wells.
- d. Add distilled water to all standard and blank wells as indicated on the label of the standard strips (A1, A2 to H1, H2).

Table 1

Table depicting an example of the arrangement of blanks, standards and samples in the microwell strips:

	1	2	3	4
A	Standard 1 (20 U/ml)	Standard 1 (20 U/ml)	Sample 1	Sample 1
B	Standard 2 (10 U/ml)	Standard 2 (10 U/ml)	Sample 2	Sample 2
C	Standard 3 (5 U/ml)	Standard 3 (5 U/ml)	Sample 3	Sample 3
D	Standard 4 (2.5 U/ml)	Standard 4 (2.5 U/ml)	Sample 4	Sample 4
E	Standard 5 (1.25 U/ml)	Standard 5 (1.25 U/ml)	Sample 5	Sample 5
F	Standard 6 (0.625 U/ml)	Standard 6 (0.625 U/ml)	Sample 6	Sample 6
G	Blank	Blank	Sample 7	Sample 7
H	Blank	Blank	Sample 8	Sample 8

- e. Add 50 µl of each prediluted **Sample**, in duplicate, to the designated wells and mix the contents.
- f. Cover with a **Plate Cover** and incubate at room temperature (18°C to 25°C) for 3 hours, if available on a microplate shaker at 200 rpm.
- g. Remove **Plate Cover** and empty wells. Wash the microwell strips 6 times with approximately 400 µl Wash Buffer per well with thorough aspiration of microwell contents between washes. Take care not to scratch the surface of the microwells.

After the last wash, tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing or place upside down on a wet absorbent paper for no longer than 15 minutes. Do not allow wells to dry.

- h. Pipette 100 µl of **TMB Substrate Solution** to all wells, including the blank wells.
- i. Incubate the microwell strips at room temperature (18° to 25 °C) for about 10 min. Avoid direct exposure to intense light.

The colour development on the plate should be monitored and the substrate reaction stopped (see point j. of this protocol) before positive wells are no longer properly recordable. Determination of the ideal time period for colour development has to be done individually for each assay.

It is recommended to add the stop solution when the highest standard has developed a dark blue colour. The colour development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.6 – 0.65.

- j. Stop the enzyme reaction by quickly pipetting 100 µl of **Stop Solution** into each well, including the blank wells. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2 - 8°C in the dark.
- k. Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the

reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the human sCD27 standards.

Note: In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.

11 Calculation of Results

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20 per cent of the mean.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the human sCD27 concentration on the abscissa. Draw a best fit curve through the points of the graph.
- To determine the concentration of circulating human sCD27 for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding human sCD27 concentration.
- ***Samples have been diluted 1:50, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 50).**
- It is suggested that each testing facility establishes a control sample of known human sCD27 concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.
- A representative standard curve is shown in Figure 4. This curve cannot be used to derive test results. Every laboratory must prepare a standard curve for each group of microwell strips assayed.

* N.B: There is a common dilution factor for samples due to the conjugate which must then be included in the calculation. The samples contribute 100 µl to the final volume per well. These 100 µl are composed of 50 µl of sample diluent plus 50 µl of the sample. This is a 1:2 dilution.

The remaining 50 µl to give 150 µl are due to the addition of 50 µl conjugate to all wells.

50 µl sample diluent and 50 µl conjugate results in 100 µl reconstitution volume, addition of 50 µl 1:25 prediluted sample (50 µl + 50 µl = 1:50 dilution)

Figure 4

Representative standard curve for human sCD27 Instant ELISA. Human sCD27 was diluted in serial 2-fold steps in sample diluent, each symbol represents the mean of 3 parallel titrations.

Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

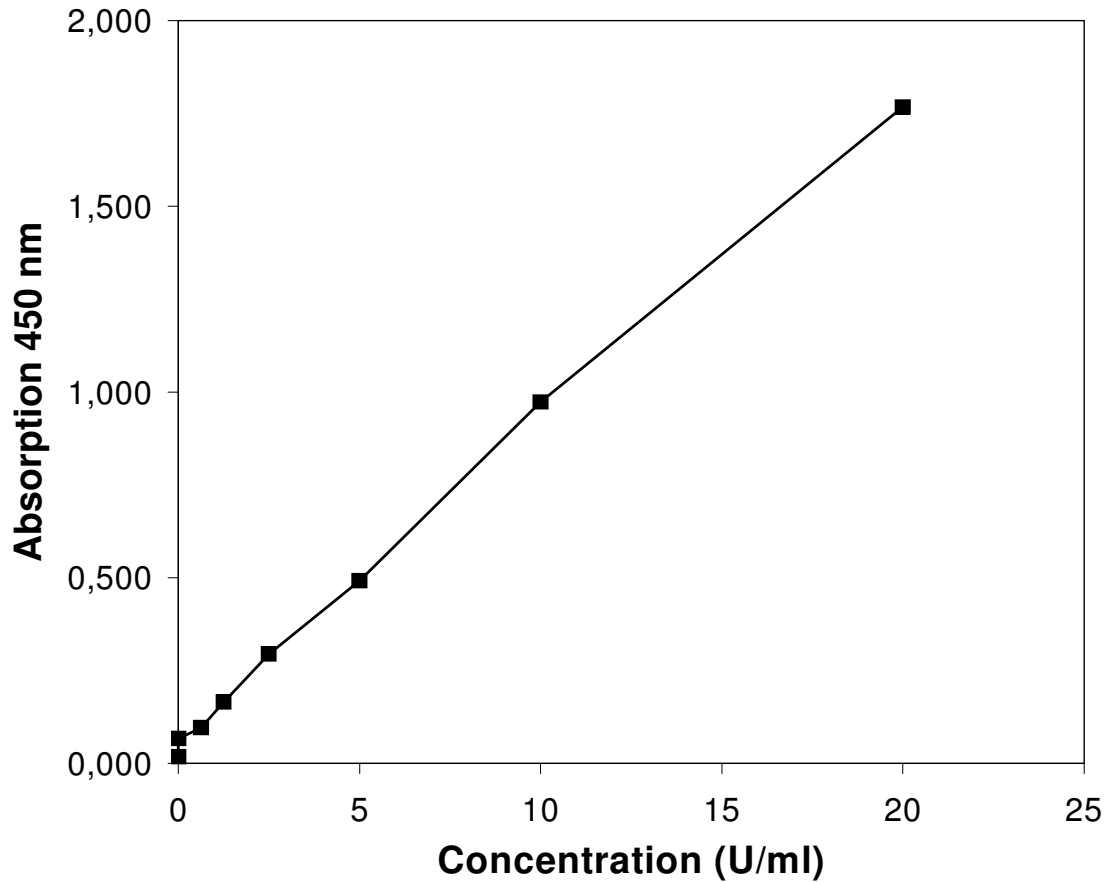


Table 2

Typical data using the human sCD27 INSTANT ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	human sCD27 Concentration (U/ml)	O.D. (450 nm)	O.D. Mean	C.V. (%)
1	20	1.805	1.767	3.0
	20	1.729		
2	10	1.022	0.973	5.1
	10	0.923		
3	5	0.533	0.492	8.4
	5	0.450		
4	2.5	0.275	0.294	6.5
	2.5	0.313		
5	1.25	0.164	0.165	0.3
	1.25	0.165		
6	0.625	0.104	0.096	8.3
	0.625	0.088		
Blank	0	0.068	0.067	
	0	0.065		

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus colour intensity. Values measured are still valid.

12 Limitations

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.
- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.
- The use of radioimmunotherapy has significantly increased the number of patients with human anti-mouse IgG antibodies (HAMA). HAMA may interfere with assays utilizing murine monoclonal antibodies leading to both false positive and false negative results. Serum samples containing antibodies to murine immunoglobulins can still be analysed in such assays when murine immunoglobulins (serum, ascitic fluid, or monoclonal antibodies of irrelevant specificity) are added to the Sample.

13 Performance Characteristics

13.1 Sensitivity

The limit of detection of human sCD27 defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be 0.2 U/ml (mean of 6 independent assays).

13.2 Reproducibility

13.2.1 Intra-assay

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 serum samples containing different concentrations of human sCD27. 2 standard curves were run on each plate. Data below show the mean human sCD27 concentration and the coefficient of variation for each sample (see Table 3). The calculated overall intra-assay coefficient of variation was 7.9%.

Table 3

The mean human sCD27 concentration and the coefficient of variation for each sample.

Positive Sample	Experiment	human sCD27 Concentration (U/ml)	Coefficient of Variation (%)
1	1	72.1	9
	2	77.6	6
	1	80.5	9
2	1	187.7	9
	2	208.8	8
	3	208.9	6
3	1	138.7	10
	2	146.8	5
	3	146.4	4
4	1	209.2	8
	2	187.0	5
	3	206.7	10
5	1	182.0	9
	2	187.5	10
	3	174.9	6
6	1	166.4	8
	2	162.4	9
	3	171.0	9
7	1	330.5	7
	2	312.9	5
	3	378.5	7
8	1	132.4	9
	2	158.5	10
	3	134.5	10

13.2.2 Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments by 3 technicians. Each assay was carried out with 6 replicates of 8 serum samples containing different concentrations of human sCD27. 2 standard curves were run on each plate. Data below (see Table 4) show the mean human sCD27 concentration and the coefficient of variation calculated on 18 determinations of each sample. The calculated overall coefficient of variation was 5.9%.

Table 4

The mean human sCD27 concentration and the coefficient of variation calculated on 18 determinations of each sample.

Sample	human sCD27 Concentration (U/ml)	Coefficient of Variation (%)
1	76.7	5.6
2	201.8	6.1
3	144.0	3.2
4	201.0	6.1
5	181.5	3.5
6	166.6	2.6
7	340.6	10.0
8	141.8	10.2

13.3 Spike Recovery

The spike recovery was evaluated by spiking 4 levels of human sCD27 into normal human serum. Recoveries were determined in 3 independent experiments with 6 replicates each. The unspiked serum was used as blank in these experiments. Average recovery ranged from 78% to 104% with an overall mean recovery of 83%.

13.4 Dilution Parallelism

4 serum samples with different levels of human sCD27 were analysed at serial 2 fold dilutions with 4 replicates each. The recovery ranged between 82.7% and 115.9% with an overall recovery of 101.7% (see Table 5).

Table 5

Sample	human sCD27 Concentration (U/ml)		% Recovery of Exp. Val.	
	Dilution	Expected Value		Observed Value
1	1:50	--	330.8	--
	1:100	165.4	191.7	115.9
	1:200	82.7	89.1	107.7
	1:400	41.4	34.2	82.7
2	1:50	--	341,0	--
	1:100	170,5	175,4	102,8%
	1:200	85,3	82,0	96,2%
	1:400	42,6	42,7	100,1%
3	1:50	--	620,6	--
	1:100	310,3	311,0	100,2%
	1:200	155,2	155,9	100,5%
	1:400	77,6	70,2	90,5%
4	1:50	--	253,5	--
	1:100	126,8	128,6	101,4%
	1:200	63,4	68,0	107,3%
	1:400	31,7	36,3	114,7%

13.5 Sample Stability

13.5.1 Freeze-Thaw Stability

Aliquots of serum samples (unspiked or spiked) were stored at -20°C and thawed 5 times, and the human sCD27 levels determined. There was no significant loss of human sCD27 immunoreactivity by freezing and thawing.

13.5.2 Storage Stability

Aliquots of serum samples (spiked or unspiked) were stored at -20 °C, 2-8 °C, room temperature (RT) and at 37 °C, and the human sCD27 level determined after 24 h. There was no significant loss of human sCD27 immunoreactivity during storage under above conditions.

13.6 Specificity

To define the specificity of this ELISA several proteins were tested for cross reactivity. There was no cross reactivity observed.

13.7 Expected Values

A panel of 16 sera from randomly selected healthy blood donors (males and females) was tested for human sCD27. The detected human sCD27 levels ranged between 269 and 424 U/ml with a mean level of 333 U/ml. The normal levels measured may vary with the sample collective used.

14 Bibliography

- 1) Agematsu K. Memory B cells and CD27. *Histol Histopathol.* 2000 Apr; 15(2): 573-6
- 2) Agematsu K, Hokiabara S, Nagumo H, Shinozaki K, Yamada S, Komiyama A. Plasma cell generation from B-lymphocytes via CS27/CD70 interaction. *Leuk Lymphoma.* 1999 Oct; 35(3-4): 219-25
- 3) Bohnhors JO, Bjorgan MB, Thoen JE, Jonsson R, Natvig JB, Thompson KM. Abnormal B cell differentiation in primary Sjogren's syndrome results in a depressed percentage of circulating memory B cells and elevated levels of soluble CD27 that correlate with Serum IgG concentration. *Clin Immunol.* 2002 Apr; 103(1): 79-88
- 4) de Jong R, Brouwer M, Kuiper HM, Hooibrink B, Miedema F, van Lier RA. Maturation- and differentiation-dependent responsiveness of human CD4+ T helper subsets. *J Immunol.* 1992 Oct 15; 149(8): 2795-802
- 5) de Jong R, Brouwer M, Hooibrink B, Van der Pouw-Kraan T, Miedema F, Van Lier RA. The CD27-subset of peripheral blood memory CD4+ lymphocytes contains functionally differentiated T lymphocytes that develop by persistent antigenic stimulation in vivo. *Eur J Immunol.* 1992 Apr ; 22(4): 993-9
- 6) Hintzen RQ, de Jong R, Lens SM, van Lier RA. CD27 : marker and mediator of T-cell activation? *Immunol Today.* 1994 Jul; 15(7): 307-11
- 7) Hintzen RQ, de Jong R, Lens SM, Brouwer M, Baars P, van Lier RA. Regulation of CD27 expression on subsets of mature T-lymphocytes. *J Immunol.* 1993 Sep 1; 151(5): 2426-35
- 8) Idriss HT, Naismith JH. TNF alpha and the TNF receptor superfamily: structure-function relationship(s). *Microsc Res Tech.* 200 Aug; 50(3): 184-95
- 9) Jacquot S. CD27/CD70 interactions regulate T dependent B cell differentiation. *Immunol Res.* 2002; 21(1): 23-30

- 10) Lens SM, de Jong R, Hintzen RQ, Koopman G, van Lier RA, van Oers RH. CD27-CD70 interaction: unravelling its implication in normal and neoplastic B-cell growth. *Leuk Lymphoma*. 1995 Jun; 18(1-2): 51-9
- 11) Lens SM, Tesselaar K, van Oers MH, van Lier RA. Control of lymphocyte function through CD27-CD70 interactions. *Semin Immunol*. 1998 Dec; 10(6): 491-9
- 12) Maurer D, Fischer GF, Fae I, Majdic O, Stuhlmeier K, Von Jeney N, Holter W, Knapp W. IgM and IgG but not cytokine secretion is restricted to the CD27+ B lymphocyte subset. *J Immunol*. 1992 Jun 15; 148(12): 3700-5
- 13) Morimoto C, Schlossman SF, P. Rambotti Lecture. Human naive and memory T cells revisited: new markers (CD31 und CD27) that help define CD4+ T cell subsets. *Clin Exp Rheumatol*. 1993 May-Jun; 11(3): 241-7
- 14) Van Baarle D, Kostense S, van Oers MH, Hamann D, Miedema F. Failing immune control as a result of impaired CD8+ T-cell maturation: CD27 might provide a clue. *Trends Immunol*. 2002 Dec; 23(12): 586-91
- 15) Wenisch C, Patruta S. Soluble CD27 in thyroid disorders. *J Lab Clin Med*. 1998 Dec; 132(6): 478-82

For literature update refer to www.bendermedsystems.com/44.html

15 Ordering Information

For orders please contact:

Europe-Headquarters

Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe
phone: +43 1 796 40 40 ext. 114
fax: +43 1 796 40 40 ext. 400
e-mail: order@bendermedsystems.com

USA

Bender MedSystems, Inc.
849 Hinckley Road
Burlingame, CA 94010, USA
toll-free phone: +1 (866) 952 2112
toll-free fax: +1 (877) 952 2112
phone: +1 (650) 952 2112
fax: +1 (650) 952 2252

For technical information please contact:

e-mail: techserv@bendermedsystems.com
www.bendermedsystems.com

Cat.No. BMS286INST human sCD27 INSTANT ELISA

16 Reagent Preparation Summary

16.1 Wash Buffer

Add **Wash Buffer Concentrate** 20 x (25 ml) to 475 ml distilled water

17 Test Protocol Summary

- Predilute sample with Sample Diluent 1:25
- Place standard strips in position A1/A2 to H1/H2.
- Add 100 µl **distilled water** to sample wells.
- Add **distilled water**, in duplicate, to all standard and blank wells as indicated on the label of the standard strips.
- Add 50 µl **Sample** to designated wells.
- Cover microwell strips and incubate 3 hours at room temperature (18° to 25°C) on a microplate shaker at 200 rpm.
- Empty and wash microwell strips 6 times with 400 µl **Wash Buffer**.
- Add 100 µl of **TMB Substrate Solution** to all wells including blank wells.
- Incubate the microwell strips for about 10 minutes at room temperature (18° to 25°C).
- Add 100 µl **Stop Solution** to all wells including blank wells.
- Blank microwell reader and measure colour intensity at 450 nm.

Note: Samples have been diluted 1:50, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 50).