

PRODUCT INFORMATION AND MANUAL

***monkey sPECAM-1
Instant ELISA***

BMS651INST

Enzyme-linked immunosorbent assay for quantitative
detection of monkey sPECAM-1.

For research use only.

Not for diagnostic or therapeutic procedures.

128 Tests

**monkey sPECAM-1
BMS651INST**



Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe

TABLE OF CONTENTS

1	Intended Use	3
2	Summary	3
3	Principles of the Test	5
4	Reagents Provided	6
5	Storage Instructions	6
6	Specimen Collection	7
7	Materials Required But Not Provided	8
8	Precautions for Use	9
9	Preparation of Reagents and Samples	11
10	Test Protocol	12
11	Calculation of Results	16
12	Limitations	19
13	Performance Characteristics	20
14	Bibliography	22
15	Ordering Information	24
16	Reagent Preparation Summary	25
17	Test Protocol Summary	26

1 Intended Use

The monkey sPECAM-1 Instant ELISA is an enzyme-linked immunosorbent assay for the quantitative detection of soluble PECAM-1 levels in cell culture supernatants, monkey serum, plasma or other body fluids. The monkey sPECAM-1 Instant ELISA is **for research use only. Not for diagnostic or therapeutic procedures.**

2 Summary

PECAM-1 (platelet endothelial cell adhesion molecule-1) also called CD31 and EndoCAM is a newly characterized adhesion molecule that belongs to the immunoglobulin superfamily [8]. PECAM-1 is a transmembrane glycoprotein with a molecular weight of approximately 130 kDa, depending on the degree of glycosylation [4].

PECAM-1 is constitutively expressed on all vascular cells and has provided a useful immunohistochemical marker of blood vessels, particularly in the setting of angiogenesis. It has also been found on platelets, monocytes, neutrophils and CD8+ T cells. Bone marrow stem cells and transformed cell lines of the myeloid and megakaryocytic lineage also express PECAM-1 [4,5]. Interestingly, PECAM-1 was also detected on solid tumor lines [9].

Recent studies suggest a role for PECAM-1 in the inflammatory process and leukocyte-endothelial interaction. The process of leukocyte emigration to the site of inflammation can be dissected into three successive stages: (1) rolling, mediated by the selectins; (2) tight adhesion mediated by ICAMs and their counter-receptors, the integrins; and (3) transmigration of leukocytes through intercellular junctions of vascular endothelial cells which requires PECAM-1 [2,7,10]. PECAM-1 appears to be able to interact both with itself (homophilic binding) and with other "non-PECAM-1" molecules (heterophilic binding).

PECAM-1 is an early and sensitive marker for tumor-induced angiogenesis [1]. Several data have suggested that PECAM-1 may be involved in the process of angiogenesis in a developing vertebrate embryo [2] as well as during metastases formation [9].

Besides the membrane-bound form of PECAM-1 a soluble form of the molecule exists, which is smaller than cell-associated PECAM-1, and contains the cytoplasmic tail. This form of soluble PECAM-1 is encoded by an alternatively spliced mRNA from which the exon containing the

transmembrane domain has been removed. Soluble PECAM-1 can be detected in plasma [6].

3 Principles of the Test

An anti-monkey sPECAM-1 monoclonal coating antibody is adsorbed onto microwells. Monkey sPECAM-1 present in the sample or standard binds to antibodies adsorbed to the microwells; an HRP-conjugated monoclonal anti-monkey sPECAM-1 antibody binds to monkey sPECAM-1 captured by the first antibody. Following incubation unbound enzyme conjugated anti-monkey sPECAM-1 is removed during a wash step and substrate solution reactive with HRP is added to the wells.

A coloured product is formed in proportion to the amount of soluble monkey sPECAM-1 present in the sample. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from seven monkey sPECAM-1 standard dilutions and monkey sPECAM-1 sample concentration determined.

Figure 1

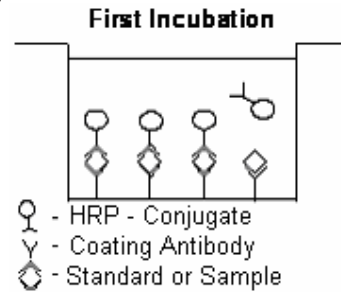


Figure 2


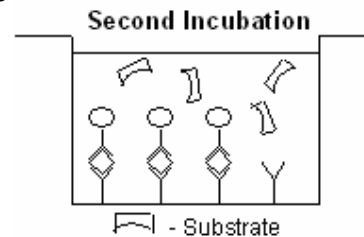


Figure 3

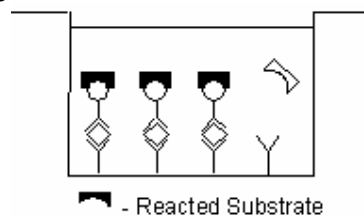

4 Reagents Provided

- 1 aluminium pouch with a **Microwell Plate coated with Monoclonal Antibody** (murine) to monkey sPECAM-1, and **HRP-Conjugate** (anti monkey sPECAM-1 monoclonal (murine) antibody), Sample Diluent lyophilized
- 2 aluminium pouches with a **monkey sPECAM-1 Standard curve (coloured)**
- 1 bottle (25 ml) **Wash Buffer Concentrate** 20x (phosphate-buffered saline with 1% Tween 20)
- 1 vial (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 1 vial (12 ml) **Sample Diluent** (Use when an external predilution of the samples is needed)
- 1 vial (15 ml) **Stop Solution** (1M Phosphoric acid)
- 2 adhesive **Plate Covers**

5 Storage Instructions

Store ELISA plate and Standard curves or whole kit at -20°C. The plate and the standard curves can also be removed, stored at -20°C, remaining kit reagents can be stored between 2° and 8°C. Expiry of the kit and reagents is stated on labels.

The expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, the reagent is not contaminated by the first handling.

6 Specimen Collection

Cell culture supernatants, monkey serum, plasma or other body fluids are suitable for use in the assay. Remove the serum or plasma from the clot or red cells as soon as possible after clotting and separation.

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples must be stored frozen at -20°C to avoid loss of bioactive monkey sPECAM-1. If samples are to be run within 24 hours, they may be stored at 2° to 8°C (for sample stability refer to 13).

Avoid repeated freeze-thaw cycles. Prior to assay, frozen serum or plasma should be brought to room temperature slowly and mixed gently.

Cell culture media without serum component are not suitable for monkey sPECAM-1 determination with the ELISA.

7 Materials Required But Not Provided

- 5 ml and 10 ml graduated pipettes
- 5 μ l to 1000 μ l adjustable single channel micropipettes with disposable tips
- 50 μ l to 300 μ l adjustable multichannel micropipette with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform linear regression analysis

8 Precautions for Use

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statements(s) for specific advice.
- Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipette by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing substrate reagent.

- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5°C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

9 Preparation of Reagents and Samples

9.1 Wash Buffer

If crystals have formed in the Wash Buffer Concentrate, warm it gently until they have completely dissolved.

Pour entire contents (25 ml) of the Wash Buffer Concentrate into a clean 500 ml graduated cylinder. Bring to final volume to 500 ml with glass-distilled or deionized water. Mix gently to avoid foaming. The pH of the final solution should adjust to 7.4.

Transfer to a clean wash bottle and store at 2° to 25°C. Please note that Wash Buffer is stable for 30 days.

10 Test Protocol

- **Use plate immediately after removal from -20°C!**
 - **Do not wait until pellets have completely dissolved before applying samples - the binding reaction in the standard strips starts immediately after addition of water!**
 - **Do not try to dissolve pellets by pipetting up and down in the wells - some parts of the pellet could stick to the tip creating high variation of results**
 - **Perform the washing step with at least 400 µl of washing buffer as stated in the manual or fill the wells completely - otherwise any pellet residues sticking to the rim of the well will not be removed and create high variation of results**
 - **Allow the washing buffer to sit in the wells for a few seconds before aspiration**
 - **Remove covers of the standard strips carefully in order that all the lyophilised pellets remain in the wells**
- a. Determine the number of microwell Strips required to test the desired number of samples plus microwell Strips for blanks and standards (coloured). Each sample, standard, blank, and optional control sample should be assayed in duplicate. Remove extra microwell Strips from holder and store in foil bag with the desiccant provided at -20°C sealed tightly. Place microwell strips containing the standard curve in position A1/A2 to H1/H2 (see Table 1).
 - b. Add 140 µl of distilled water to the sample wells.
 - c. Add distilled water to all standard and blank wells as indicated on the label of the standard strips (A1, A2 to H1, H2).

Table 1

Table depicting an example of the arrangement of blanks, standards and samples in the microwell strips:

	1	2	3	4
A	Standard 1 (30 U/ml)	Standard 1 (30 U/ml)	Sample 1	Sample 1
B	Standard 2 (15 U/ml)	Standard 2 (15 U/ml)	Sample 2	Sample 2
C	Standard 3 (7.5 U/ml)	Standard 3 (7.5 U/ml)	Sample 3	Sample 3
D	Standard 4 (3.75 U/ml)	Standard 4 (3.75 U/ml)	Sample 4	Sample 4
E	Standard 5 (1.9 U/ml)	Standard 5 (1.9 U/ml)	Sample 5	Sample 5
F	Standard 6 (0.95 U/ml)	Standard 6 (0.95 U/ml)	Sample 6	Sample 6
G	Standard 7 (0.48 U/ml)	Standard 7 (0.48 U/ml)	Sample 7	Sample 7
H	Blank	Blank	Sample 8	Sample 8

- d. Add 10 µl of each **Sample**, in duplicate, to the designated wells and mix the contents.
- e. Cover with a **Plate Cover** and incubate at room temperature (18°C to 25°C) for 4 hours, if available on a microplate shaker at 100 rpm.
- f. Remove **Plate Cover** and empty wells. Wash the microwell strips 3 times with approximately 400 µl Wash Buffer per well with thorough aspiration of microwell contents between washes. Take care not to scratch the surface of the microwells.

After the last wash, tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing or place upside down on a wet absorbent paper for no longer than 15 minutes. Do not allow wells to dry.

- g. Pipette 100 µl of **TMB Substrate Solution** to all wells, including the blank wells.
- h. Incubate the microwell strips at room temperature (18° to 25 °C) for about 10 min. Avoid direct exposure to intense light.

The colour development on the plate should be monitored and the substrate reaction stopped (see point i. of this protocol) before positive wells are no longer properly recordable. Determination of the ideal time period for colour development has to be done individually for each assay.

It is recommended to add the stop solution when the highest standard has developed a dark blue colour. The colour development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.6 – 0.65.

- i. Stop the enzyme reaction by quickly pipetting 100 µl of **Stop Solution** into each well, including the blank wells. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2 - 8°C in the dark.
- j. Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the

reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the monkey sPECAM-1 standards.

Note: In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.

11 Calculation of Results

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20 per cent of the mean.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the monkey sPECAM-1 concentration on the abscissa. Draw a best fit curve through the points of the graph.
- To determine the concentration of circulating monkey sPECAM-1 for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding monkey sPECAM-1 concentration.
- ***Samples have been diluted 1:10, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 10).**
- It is suggested that each testing facility establishes a control sample of known monkey sPECAM-1 concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.
- A representative standard curve is shown in Figure 4. This curve cannot be used to derive test results. Every laboratory must prepare a standard curve for each group of microwell strips assayed.

* N.B: There is a common dilution factor for samples due to the conjugate which must then be included in the calculation. The samples contribute 100 µl to the final volume per well. These 100 µl are composed of 90 µl of sample diluent plus 10 µl of the sample. This is a 1:10 dilution.

The remaining 50 µl to give 150 µl are due to the addition of 50 µl conjugate to all wells.

90 µl sample diluent and 50 µl conjugate results in 140 µl reconstitution volume, addition of 10 µl sample (90 µl + 10 µl = 1:10 dilution)

Figure 4

Representative standard curve for monkey sPECAM-1 Instant ELISA. Monkey sPECAM-1 was diluted in serial 2-fold steps in Sample Diluent, each symbol represents the mean of 3 parallel titrations. Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

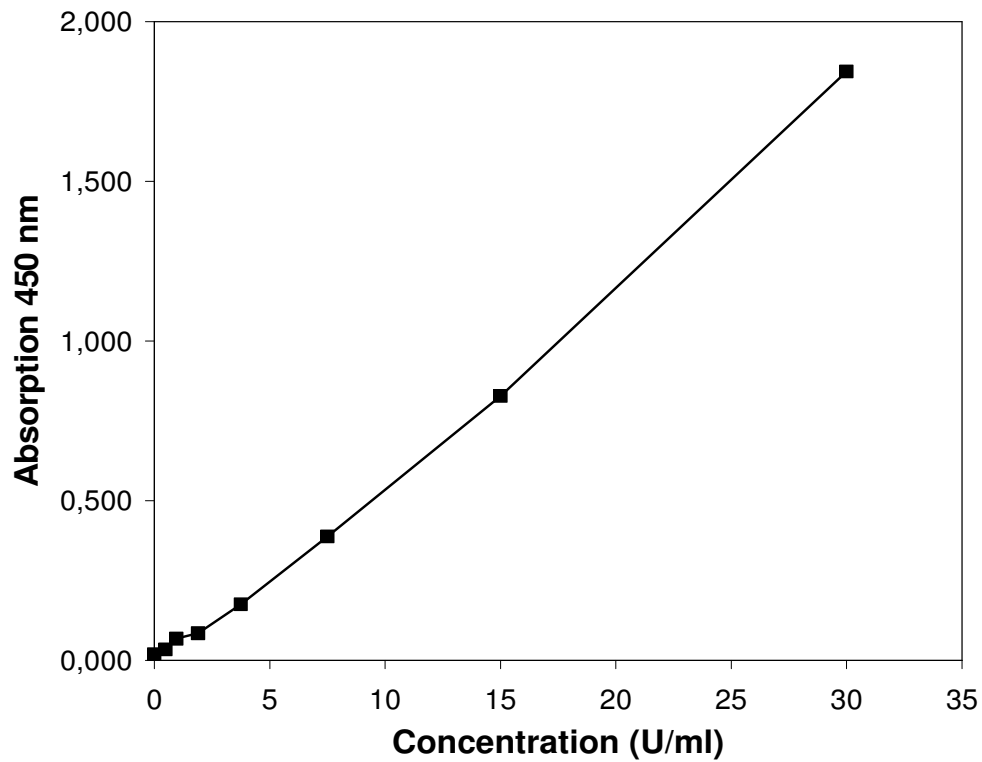


Table 2

Typical data using the monkey sPECAM-1 INSTANT ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	monkey sPECAM-1 Concentration (U/ml)	O.D. (450 nm)	O.D. Mean
1	30	1.852	1.844
	30	1.835	
2	15	0.830	0.828
	15	0.825	
3	7.5	0.381	0.388
	7.5	0.395	
4	3.75	0.180	0.176
	3.75	0.172	
5	1.9	0.081	0.085
	1.9	0.088	
6	0.95	0.067	0.068
	0.95	0.069	
7	0.48	0.037	0.034
	0.48	0.031	
Blank	0	0.019	0.019
	0	0.018	

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus colour intensity. Values measured are still valid.

12 Limitations

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.
- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.
- The use of mouse IgG antibodies or mouse IgG antibody coupled molecules in the course of animal experiments has significantly increased the number of animals with anti-mouse IgG antibodies (equivalent to human anti-mouse IgG antibodies - HAMA). These anti-mouse IgG antibodies may interfere with assays utilizing murine monoclonal antibodies leading to both false positive and false negative results. Serum samples containing antibodies to murine immunoglobulins can still be analysed in such assays when murine immunoglobulins (serum, ascitic fluid, or monoclonal antibodies of irrelevant specificity) are added to the Sample.

13 Performance Characteristics

13.1 Sensitivity

The limit of detection of monkey sPECAM-1 defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be 0.06 U/ml (mean of 10 independent assays).

13.2 Reproducibility

13.2.1 Intra-assay

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 samples containing different concentrations of monkey sPECAM-1. 2 standard curves were run on each plate. The calculated overall intra-assay coefficient of variation was <10%.

13.2.2 Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments by 3 technicians. Each assay was carried out with 6 replicates of 8 samples containing different concentrations of monkey sPECAM-1. 2 standard curves were run on each plate. The calculated overall inter-assay coefficient of variation was 10%.

13.3 Spike Recovery

The spike recovery was evaluated by spiking 4 levels of monkey sPECAM-1 into samples (serum matrix). Recoveries were determined in 3 independent experiments with 6 replicates each. The unspiked sample was used as blank in these experiments. The overall mean recovery was 105%.

13.4 Dilution Parallelism

4 samples with different levels of monkey sPECAM-1 were analysed at serial 2 fold dilutions with 4 replicates each. The overall mean recovery was 90%.

13.5 Sample Stability

13.5.1 Freeze-Thaw Stability

Aliquots of samples (unspiked or spiked) were stored at -20 °C and thawed 5 times, and the monkey sPECAM-1 levels determined. There was no significant loss of monkey sPECAM-1 immunoreactivity detected by freezing and thawing.

13.5.2 Storage Stability

Aliquots of samples (spiked or unspiked) were stored at -20 °C, 2-8 °C, room temperature (RT) and at 37 °C, and the monkey sPECAM-1 level determined after 24, 48 and 96 h. There was no significant loss of monkey sPECAM-1 immunoreactivity detected during storage under above conditions.

14 Bibliography

- 1) Berger, R., Albelda, S.M., Berd, D., Ioffreda, M., Whitaker, D. and Murphy, G.F. (1993). Expression of platelet-endothelial cell adhesion molecule-1 (PECAM-1) during melanoma-induced angiogenesis in vivo. *J. Cutan. Pathol.* 20, 399-406.
- 2) Bogen, S., Pak, J., Garifallou, M., Deng, X. and Muller, W.A. (1994). Monoclonal antibody to murine PECAM-1 (CD31) blocks acute inflammation in vivo. *J. Exp. Med.* 179, 1059-1064.
- 3) Buck, C.A., Baldwin, H.S., DeLisser, H., Mickanin, C., Shen, H.M., Kennedy, G., Chen, A., Edelman, J.M. and Albelda, S.M. (1993). Cell adhesion receptors and early mammalian heart development: an overview. *C. R. Acad. Sci. III.* 316, 838-859.
- 4) DeLisser, H.M., Newman, P.J. and Albelda, S.M. (1993). Platelet endothelial cell adhesion molecule (CD31). *Curr. Top. Microbiol. Immunol.* 184, 37-45.
- 5) DeLisser, H.M., Newman, P.J. and Albelda, S.M. (1994). Molecular and functional aspects of PECAM-1/CD31. *Immunol. Today* 15,
- 6) Goldberger, A., Middleton, K.A., Oliver, J.A., Paddock, C., Yan, H.C., DeLisser, H.M., Albelda, S.M. and Newman, P.J. (1994). Biosynthesis and processing of the cell adhesion molecule PECAM-1 includes production of a soluble form. *J. Biol. Chem.* 269, 17183-17191.
- 7) Muller, W.A., Weigl, S.A., Deng, X. and Phillips, D.M. (1993). PECAM-1 is required for transendothelial migration of leukocytes. *J. Exp. Med.* 178, 449-460.
- 8) Newman, P.J., Berndt, M.C., Gorski, J., White II, G.C., Lyman, S., Paddock, C. and Muller, W.A. (1990). PECAM-1 (CD31) Cloning and Relation to Adhesion Molecules of the Immunoglobulin Gene Superfamily. *Science* 247, 1219-1222.
- 9) Tang, D.G., Chen, Y.Q., Newman, P.J., Shi, L., Gao, X., Diglio, C.A. and Honn, K.V. (1993). Identification of PECAM-1 in solid tumor cells and its potential involvement in tumor cell adhesion to endothelium. *J. Biol. Chem.* 268, 22883-22894.

- 10) Vaporciyan, A.A., DeLisser, H.M., Yan, H.C., Mendiguren, I.I., Thom, S.R., Jones, M.L., Ward, P.A. and Albelda, S.M. (1993). Involvement of platelet-endothelial cell adhesion molecule-1 in neutrophil recruitment in vivo. *Science* 262, 1580-1582.

For literature update refer to **www.bendermedsystems.com/44.html**

15 Ordering Information

For orders please contact:

Europe-Headquarters

Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe
phone: +43 1 796 40 40 ext. 114
fax: +43 1 796 40 40 ext. 400
e-mail: order@bendermedsystems.com

USA

Bender MedSystems, Inc.
849 Hinckley Road
Burlingame, CA 94010, USA
toll-free phone: +1 (866) 952 2112
toll-free fax: +1 (877) 952 2112
phone: +1 (650) 952 2112
fax: +1 (650) 952 2252

For technical information please contact:

e-mail: techserv@bendermedsystems.com
www.bendermedsystems.com

Cat.No. BMS651INST monkey sPECAM-1 INSTANT ELISA

16 Reagent Preparation Summary

16.1 Wash Buffer

Add **Wash Buffer Concentrate** 20 x (25 ml) to 475 ml distilled water

17 Test Protocol Summary

- Place standard strips in position A1/A2 to H1/H2.
- Add 140 µl **distilled water** to sample wells.
- Add **distilled water**, in duplicate, to all standard and blank wells as indicated on the label of the standard strips.
- Add 10 µl **Sample** to designated wells.
- Cover microwell strips and incubate 4 hours at room temperature (18° to 25°C) on a microplate shaker at 100 rpm.
- Empty and wash microwell strips 3 times with 400 µl **Wash Buffer**.
- Add 100 µl of **TMB Substrate Solution** to all wells including blank wells.
- Incubate the microwell strips for about 10 minutes at room temperature (18° to 25°C).
- Add 100 µl **Stop Solution** to all wells including blank wells.
- Blank microwell reader and measure colour intensity at 450 nm.

Note: Samples have been diluted 1:10, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 10).