

PRODUCT INFORMATION AND MANUAL

***mouse RANTES
Instant ELISA***

BMS6009INST

Enzyme-linked immunosorbent assay for quantitative
detection of mouse RANTES.

For research use only.

Not for diagnostic or therapeutic procedures.

128 Tests

**mouse RANTES
BMS6009INST**



Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe

TABLE OF CONTENTS

1	Intended Use	3
2	Summary	3
3	Principles of the Test	4
4	Reagents Provided	5
5	Storage Instructions	5
6	Specimen Collection	6
7	Materials Required But Not Provided	7
8	Precautions for Use	8
9	Preparation of Reagents and Samples	10
10	Test Protocol	11
11	Calculation of Results	15
12	Limitations	18
13	Performance Characteristics	19
14	Bibliography	23
15	Ordering Information	26
16	Reagent Preparation Summary	27
17	Test Protocol Summary	28

1 Intended Use

The mouse RANTES Instant ELISA is an enzyme-linked immunosorbent assay for the quantitative detection of murine RANTES levels in cell culture supernatants, murine serum and plasma. The mouse RANTES Instant ELISA is **for research use only. Not for diagnostic or therapeutic procedures.**

2 Summary

RANTES (Regulated upon Activation, Normal T cell Expressed and presumably Secreted), also known as CCL5, is a member of the CC chemokine family of inflammatory and immunoregulatory cytokines.

Similar to other CC chemokines RANTES is a potent chemoattractant for a number of different cell types including monocytes (9), unstimulated CD4⁺/CD45RO⁺ memory T-cells, stimulated CD4⁺ and CD8⁺ T cells with naïve and memory phenotypes (21, 23), NK cells (13), eosinophils (10), basophils (4), dendritic cells (12) and mast cells (16). The activity of RANTES however is not restricted merely to chemotaxis, moreover it is a powerful activator for T-cells, monocytes, neutrophils and NK cells, a feature potentially relevant in a range of inflammatory disorders (1, 2, 3).

Several CC chemokine receptors, including CCR-1, CCR-3, CCR-4 and CCR-5, have been shown to bind RANTES and subsequently to transduce a signal by increasing the intracellular calcium ion level (6, 17, 18, 19).

Mouse RANTES cDNA encodes a 91 amino acid residue precursor protein that is cleaved to generate a 68 aa residue mature protein (22). At the amino acid sequence level, mouse RANTES is 84% identical to human RANTES.

RANTES, amongst other chemokines, has been shown to play an important role in atherosclerosis (5), acute and chronic allergic inflammation (7, 14, 20), HIV suppression (8), food allergy (11) and macrophage accumulation in the uterus during pregnancy (24). RANTES deficient mice showed impaired T-cell function (15)

3 Principles of the Test

An anti-mouse RANTES polyclonal coating antibody is adsorbed onto microwells. Mouse RANTES present in the sample or standard binds to antibodies adsorbed to the microwells; a biotin-conjugated polyclonal anti-mouse RANTES antibody binds to mouse RANTES captured by the first antibody. Streptavidin-HRP binds to the biotin conjugated anti-mouse RANTES.

Following incubation unbound biotin conjugated anti mouse RANTES and Streptavidin-HRP is removed during a wash step, and substrate solution reactive with HRP is added to the wells.

A coloured product is formed in proportion to the amount of soluble mouse RANTES present in the sample. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from seven mouse RANTES standard dilutions and mouse RANTES sample concentration determined.

Figure 1

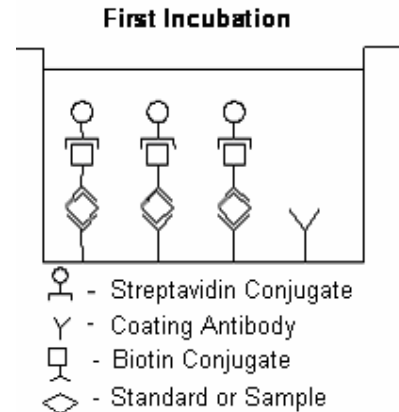


Figure 2


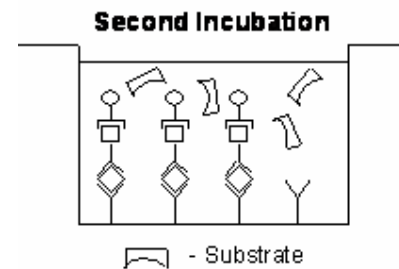


Figure 3

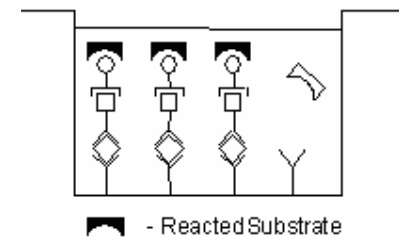

4 Reagents Provided

- 1 aluminium pouch with a **Microwell Plate coated with Polyclonal Antibody** to mRANTES, **Biotin-Conjugate** (anti-mRANTES polyclonal antibody), **Streptavidin-HRP** and Sample diluent, lyophilized
- 2 aluminium pouches with a **mouse RANTES Standard curve (coloured)**
- 1 bottle (25 ml) **Wash Buffer Concentrate** 20x (phosphate-buffered saline with 1% Tween 20)
- 1 vial (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 1 vial (12 ml) **Sample Diluent** (Use when an external predilution of the samples is needed)
- 1 vial (15 ml) **Stop Solution** (1M Phosphoric acid)
- 2 adhesive **Plate Covers**

5 Storage Instructions

Store ELISA plate and Standard curves or whole kit at -20 °C. The plate and the standard curves can also be removed, stored at -20 °C, remaining kit reagents can be stored between 2 ° and 8 °C. Expiry of the kit and reagents is stated on labels.

The expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, the reagent is not contaminated by the first handling.

6 Specimen Collection

Cell culture supernatants, murine serum, EDTA or heparinized plasma are suitable for use in the assay. Remove the serum or plasma from the clot or red cells as soon as possible after clotting and separation.

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples must be stored frozen at -20°C to avoid loss of bioactive mouse RANTES. If samples are to be run within 24 hours, they may be stored at 2° to 8°C (for sample stability refer to 13).

Avoid repeated freeze-thaw cycles. Prior to assay, frozen serum or plasma should be brought to room temperature slowly and mixed gently.

7 Materials Required But Not Provided

- 5 ml and 10 ml graduated pipettes
- 5 μ l to 1000 μ l adjustable single channel micropipettes with disposable tips
- 50 μ l to 300 μ l adjustable multichannel micropipette with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform linear regression analysis

8 Precautions for Use

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statements(s) for specific advice.
- Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipette by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing substrate reagent.

- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5 °C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

9 Preparation of Reagents and Samples

9.1 Wash Buffer

If crystals have formed in the Wash Buffer Concentrate, warm it gently until they have completely dissolved.

Pour entire contents (25 ml) of the Wash Buffer Concentrate into a clean 500 ml graduated cylinder. Bring to final volume to 500 ml with glass-distilled or deionized water. Mix gently to avoid foaming. The pH of the final solution should adjust to 7.4.

Transfer to a clean wash bottle and store at 2° to 25°C. Please note that Wash Buffer is stable for 30 days.

10 Test Protocol

- **Use plate immediately after removal from -20°C!**
 - **Do not wait until pellets have completely dissolved before applying samples - the binding reaction in the standard strips starts immediately after addition of water!**
 - **Do not try to dissolve pellets by pipetting up and down in the wells - some parts of the pellet could stick to the tip creating high variation of results**
 - **Perform the washing step with at least 400 µl of washing buffer as stated in the manual or fill the wells completely - otherwise any pellet residues sticking to the rim of the well will not be removed and create high variation of results**
 - **Allow the washing buffer to sit in the wells for a few seconds before aspiration**
 - **Remove covers of the standard strips carefully in order that all the lyophilised pellets remain in the wells**
- a. Determine the number of microwell Strips required to test the desired number of samples plus microwell Strips for blanks and standards (coloured). Each sample, standard, blank, and optional control sample should be assayed in duplicate. Remove extra microwell Strips from holder and store in foil bag with the desiccant provided at -20°C sealed tightly. Place microwell strips containing the standard curve in position A1/A2 to H1/H2 (see Table 1).
 - b. Add 100 µl of distilled water to the sample wells.
 - c. Add distilled water to all standard and blank wells as indicated on the label of the standard strips (A1, A2 to H1, H2).

Table 1

Table depicting an example of the arrangement of blanks, standards and samples in the microwell strips:

	1	2	3	4
A	Standard 1 (500 pg/ml)	Standard 1 (500 pg/ml)	Sample 1	Sample 1
B	Standard 2 (250 pg/ml)	Standard 2 (250 pg/ml)	Sample 2	Sample 2
C	Standard 3 (125 pg/ml)	Standard 3 (125 pg/ml)	Sample 3	Sample 3
D	Standard 4 (62.5 pg/ml)	Standard 4 (62.5 pg/ml)	Sample 4	Sample 4
E	Standard 5 (31.3 pg/ml)	Standard 5 (31.3 pg/ml)	Sample 5	Sample 5
F	Standard 6 (15.6 pg/ml)	Standard 6 (15.6 pg/ml)	Sample 6	Sample 6
G	Standard 7 (7.8 pg/ml)	Standard 7 (7.8 pg/ml)	Sample 7	Sample 7
H	Blank	Blank	Sample 8	Sample 8

- d. Add 50 µl of each **Sample**, in duplicate, to the designated wells and mix the contents.
- e. Cover with a **Plate Cover** and incubate at room temperature (18°C to 25°C) for 3 hours, if available on a microplate shaker at 200 rpm.
- f. Remove **Plate Cover** and empty wells. Wash the microwell strips 6 times with approximately 400 µl Wash Buffer per well with thorough aspiration of microwell contents between washes. Take care not to scratch the surface of the microwells.

After the last wash, tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing or place upside down on a wet absorbent paper for no longer than 15 minutes. Do not allow wells to dry.

- g. Pipette 100 µl of **TMB Substrate Solution** to all wells, including the blank wells.
- h. Incubate the microwell strips at room temperature (18° to 25°C) for about 10-20 min. Avoid direct exposure to intense light.

The colour development on the plate should be monitored and the substrate reaction stopped (see point i. of this protocol) before positive wells are no longer properly recordable. Determination of the ideal time period for colour development has to be done individually for each assay.

It is recommended to add the stop solution when the highest standard has developed a dark blue colour. The colour development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.6 – 0.65.

- i. Stop the enzyme reaction by quickly pipetting 100 µl of **Stop Solution** into each well, including the blank wells. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2 - 8°C in the dark.
- j. Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the

reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the mouse RANTES standards.

Note: In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.

11 Calculation of Results

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20 per cent of the mean.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the mouse RANTES concentration on the abscissa. Draw a best fit curve through the points of the graph.
- To determine the concentration of circulating mouse RANTES for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding mouse RANTES concentration.
- ***Samples have been diluted 1:2, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 2).**
- It is suggested that each testing facility establishes a control sample of known mouse RANTES concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.
- A representative standard curve is shown in Figure 4. This curve cannot be used to derive test results. Every laboratory must prepare a standard curve for each group of microwell strips assayed.

* N.B: There is a common dilution factor for samples due to the conjugate which must then be included in the calculation. The samples contribute 100 µl to the final volume per well. These 100 µl are composed of 50 µl of sample diluent plus 50 µl of the sample. This is a 1:2 dilution.

The remaining 50 µl to give 150 µl are due to the addition of 50 µl conjugate to all wells.

50 µl sample diluent and 50 µl conjugate results in 100 µl reconstitution volume, addition of 50 µl sample (50 µl + 50 µl = 1:2 dilution)

Figure 4

Representative standard curve for mouse RANTES Instant ELISA. Mouse RANTES was diluted in serial 2-fold steps in Sample Diluent, each symbol represents the mean of 3 parallel titrations. Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

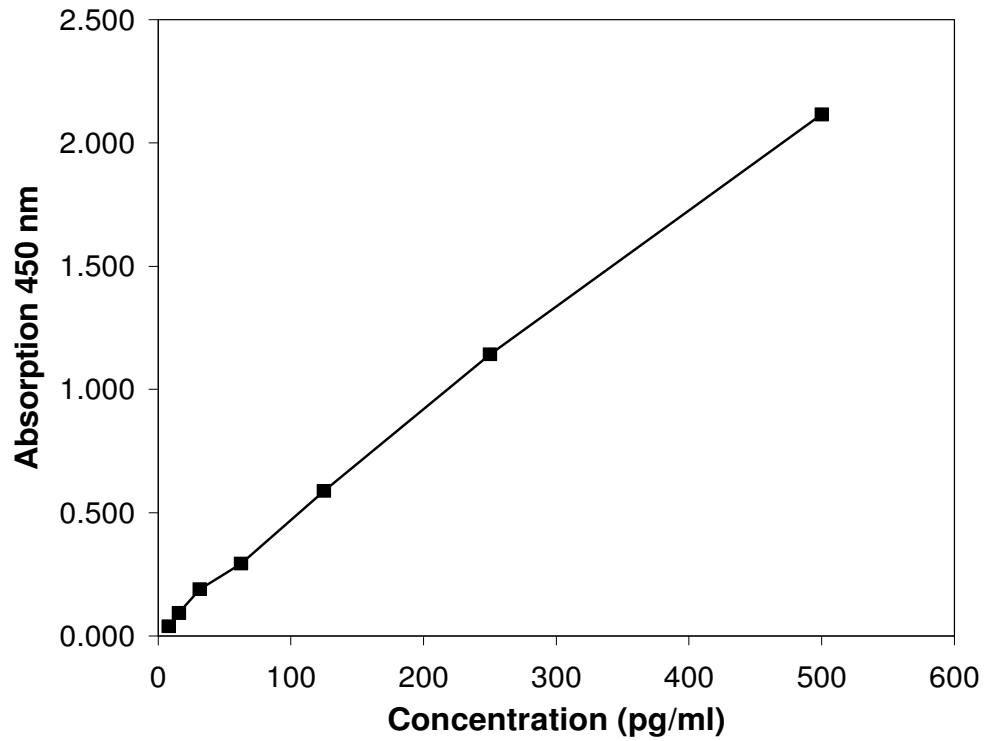


Table 2

Typical data using the mouse RANTES INSTANT ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	mouse RANTES Concentration (pg/ml)	O.D. (450 nm)	O.D. Mean	C.V. (%)
1	500	2.252	2.160	6.1
	500	2.067		
2	250	1.271	1.187	10.0
	250	1.103		
3	125	0.605	0.632	6.0
	125	0.659		
4	62.5	0.331	0.337	2.3
	62.5	0.342		
5	31.3	0.241	0.234	4.5
	31.3	0.226		
6	15.6	0.136	0.137	0.5
	15.6	0.137		
7	7.8	0.080	0.083	4.3
	7.8	0.085		
Blank	0	0.042	0.044	4.9
	0	0.045		

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus colour intensity. Values measured are still valid.

12 Limitations

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.
- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.

13 Performance Characteristics

13.1 Sensitivity

The limit of detection of mouse RANTES defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be 1.2 pg/ml (mean of 6 independent assays).

13.2 Reproducibility

13.2.1 Intra-assay

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 4 replicates of 8 samples containing different concentrations of mouse RANTES. 2 standard curves were run on each plate. Data below show the mean mouse RANTES concentration and the coefficient of variation for each sample (see Table 3). The calculated overall intra-assay coefficient of variation was 6.2%.

Table 3

The mean mouse RANTES concentration and the coefficient of variation for each sample.

Positive Sample	Experiment	mouse RANTES Concentration (pg/ml)	Coefficient of Variation (%)
1	1	820	5
	2	855	5
	3	753	7
2	1	456	2
	2	552	9
	3	387	6
3	1	266	4
	2	268	13
	3	189	6
4	1	144	9
	2	156	6
	3	107	11
5	1	1276	1
	2	1210	9
	3	1340	6
6	1	627	4
	2	548	4
	3	678	5
7	1	312	6
	2	295	4
	3	262	9
8	1	157	5
	2	148	2
	3	127	10

13.2.2 Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments by 3 technicians. Each assay was carried out with 4 replicates of 8 samples containing different concentrations of mouse RANTES. 2 standard curves were run on each plate. Data below (see Table 4) show the mean mouse RANTES concentration and the coefficient of variation calculated on 12 determinations of each sample. The calculated overall coefficient of variation was 12.1%.

Table 4

The mean mouse RANTES concentration and the coefficient of variation calculated on 18 determinations of each sample.

Sample	mouse RANTES Concentration (pg/ml)	Coefficient of Variation (%)
1	809	6.4
2	465	17.8
3	241	18.8
4	135	18.9
5	1275	5.1
6	618	10.6
7	290	8.7
8	144	10.5

13.3 Spike Recovery

The spike recovery was evaluated by spiking 4 levels of recombinant mouse RANTES into mouse serum. Recoveries were determined in 3 independent experiments with 2 replicates each. The unspiked serum was used as blank in these experiments. The overall mean recovery was 77%.

13.4 Dilution Parallelism

4 serum samples with different levels of mouse RANTES were analysed at serial 2 fold dilutions with 4 replicates each. The recovery ranged between 87% and 116% with an overall recovery of 102%.

13.5 Sample Stability

13.5.1 Freeze-Thaw Stability

Aliquots of serum samples (unspiked or spiked) were stored at -20 °C and thawed 5 times, and the mouse RANTES levels determined. There was a **significant increase** of mouse RANTES immunoreactivity detected. Therefore samples should be stored at -20 °C in aliquots and **thawed only once**.

13.5.2 Storage Stability

Aliquots of serum samples (spiked or unspiked) were stored at -20 °C, 2-8 °C, room temperature (RT) and at 37 °C, and the mouse RANTES level determined after 24 h. There was no significant loss of mouse RANTES immunoreactivity detected. during storage under above conditions.

13.6 Specificity

The assay detects both natural and recombinant mouse RANTES. To define the specificity of this ELISA several proteins were tested for cross reactivity. There was no cross reactivity observed.

13.7 Expected Values

A panel of sera from randomly selected mice was tested for mouse RANTES. The detected mouse RANTES levels ranged between 34 and 118 pg/ml with an overall mean level of 73 pg/ml.

14 Bibliography

- 1) Appay V, Rowland-Jones SL. RANTES: a versatile and controversial chemokine. *Trends Immunol.* 2001 Feb; 22(2): 83-7.
- 2) Appay V, Brown A, Cribbes S, Randle E, Czaplewski LG. Aggregation of RANTES is responsible for its inflammatory properties. Characterization of nonaggregating, noninflammatory RANTES mutants. *J Biol Chem.* 1999 Sept; 274(39): 217505-12.
- 3) Appay V, Dunbar PR, Cerundolo V, McMichael A, Czaplewski L, Rowland-Jones S. RANTES activates antigen-specific cytotoxic T lymphocytes in a mitogen-like manner through cell surface aggregation. *Int Immunol.* 2000 Aug; 12(8): 1173-82.
- 4) Bischoff SC, Krieger M, Brunner T, Rot A, von Tscharner V, Baggiolini M, Dahinden CA. RANTES and related chemokines activate human basophil granulocytes through different G protein-coupled receptors. *Eur J Immunol.* 1993 Mar; 23(3): 761-7.
- 5) Bursill CA, Channon KM, Greaves DR. The role of chemokines in atherosclerosis: recent evidence from experimental models and population genetics. *Curr Opin Lipidol.* 2004 Apr; 15(2): 145-9.
- 6) Combadiere C, Ahuja SK, Tiffany HL, Murphy PM. Cloning and functional expression of CC CKR5, a human monocyte CC chemokine receptor selective for MIP-1(alpha), MIP-1(beta) and RANTES. *J Leukoc Biol.* 1996 Jul; 60(1): 147-52.
- 7) Conti P, DiGioacchino M. MCP-1 and RANTES are mediators of acute and chronic inflammation. *Allergy Asthma Proc.* 2001 May-Jun; 22(3) : 133-7.
- 8) DeVico A, Gallo RC. Control of HIV-1 infection by soluble factors of the immune response. *Nature Reviews Microbiology* 2, 401-413 (2004).
- 9) Fine JS, Byrnes HD, Zavodny PJ, Hipkin RW. Evaluation of signal transduction pathways in chemoattractant-induced human monocyte chemotaxis. *Inflammation.* 2001 Apr; 25(2): 61-7.

- 10) Kameyoshi Y, Dorschner A, Mallet AI, Christophers E, Schroder JM. Cytokine RANTES released by thrombin-stimulated platelets is a potent attractant for human eosinophils. *J Exp Med*. 1992 Aug; 176(2): 587-92.
- 11) Lee JB, Matsumoto T, Shin YO, Yang HM, Min YK, Timothy O, Bae SJ, Quan FS. The role of RANTES in a murine model of food allergy. *Immunol Invest*. 2004 Feb; 33(1): 27-38.
- 12) Lin CL, Suri RM, Rahdon RA, Austyn JM, Roake JA. Dendritic cell chemotaxis and transendothelial migration are induced by distinct chemokines and are regulated on maturation. *Eur J Immunol*. 1998 Dec ; 28(12) : 4114-22.
- 13) Loetscher P, Seitz M, Clark-Lewis I, Baggiolini M, Moser B. Activation of NK cells by CC chemokines. Chemotaxis, Ca²⁺ mobilization, and enzyme release. *J Immunol*. 1996 Jan; 156(1): 322-7.
- 14) Lukacs NW, Standiford TJ, Chensue SW, Kunkel RG, Strieter RM, Kunkel SL. C-C chemokine-induced eosinophil chemotaxis during allergic airway inflammation. *J Leukoc Biol*. 1996 Nov; 60(5): 573-8.
- 15) Makino Y, Cook DN, Smithies O, Hwang OY, Neilson EG, Turka LA, Sato H, Wells AD, Danoff TM. Impaired T cell function in RANTES-deficient mice. *Clin Immunol*. 2002 Mar; 102(3): 302-9.
- 16) Mattoli S, Ackerman V, Vittori E, Marini M. Mast cell chemotactic activity of RANTES. *Biochem Biophys Res Commun*. 1995 Apr; 209(1): 316-21.
- 17) Neote K, Digregorio D, Mak JY, Horuk R, Schall TJ. Molecular cloning, functional expression, and signaling characteristics of a C-C chemokine receptor. *Cell*. 1993 Feb; 72(3): 415-25.
- 18) Ponath PD, Qin S, Post TW, Wang J, Wu L, Gerard NP, Newman W, Gerard C, Mackay CR. Molecular cloning and characterization of a human eotaxin receptor expressed selectively on eosinophils. *J Exp Med*. 1996 Jun; 183(6): 2437-48.
- 19) Power CA, Meyer A, Nemeth K, Bacon KB, Hoogewerf AJ, Proudfoot AE, Wells TN. Molecular cloning and functional expression of a novel CC chemokine receptor cDNA from a human basophilic cell line. *J Biol Chem*. 1995 Aug; 270(33): 19495-500.

- 20) Romagnani S. Cytokines and chemoattractants in allergic inflammation. *Mol Immunol*. 2002 May; 38(12-13): 881-5.
- 21) Schall TJ, Bacon K, Toy KJ, Goeddel DV. Selective attraction of monocytes and T lymphocytes of the memory phenotype by cytokine RANTES. *Nature*. 1990 Oct; 347(6294): 669-71.
- 22) Schall TJ, Jongstra J, Dyer BJ, Jorgensen J, Clayberger C, Davis MM, Krensky AM. A human T cell-specific molecule is a member of a new gene family. *J Immunol*. 1988 Aug; 141(3): 1018-25.
- 23) Taub DD, Oppenheim JJ. Chemokines, inflammation and the immune system. *Ther Immunol*. 1994 Aug; 1(4): 229-46.
- 24) Wood GW, Hausmann EH, Kanakaraj K. Expression and regulation of chemokine genes in the mouse uterus during pregnancy. *Cytokine*. 1999 Dec; 11(12): 1038-45.

For literature update refer to **www.bendermedsystems.com/44.html**

15 Ordering Information

For orders please contact:

Europe-Headquarters

Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe
phone: +43 1 796 40 40 ext. 114
fax: +43 1 796 40 40 ext. 400
e-mail: order@bendermedsystems.com

USA

Bender MedSystems, Inc.
849 Hinckley Road
Burlingame, CA 94010, USA
toll-free phone: +1 (866) 952 2112
toll-free fax: +1 (877) 952 2112
phone: +1 (650) 952 2112
fax: +1 (650) 952 2252

For technical information please contact:

e-mail: techserv@bendermedsystems.com
www.bendermedsystems.com

Cat.No. BMS6009INST mouse RANTES INSTANT ELISA

16 Reagent Preparation Summary

16.1 Wash Buffer

Add **Wash Buffer Concentrate** 20 x (25 ml) to 475 ml distilled water

17 Test Protocol Summary

- Place standard strips in position A1/A2 to H1/H2.
- Add 100 µl **distilled water** to sample wells.
- Add **distilled water**, in duplicate, to all standard and blank wells as indicated on the label of the standard strips.
- Add 50 µl **Sample** to designated wells.
- Cover microwell strips and incubate 3 hours at room temperature (18° to 25°C) on a microplate shaker at 200 rpm.
- Empty and wash microwell strips 6 times with 400 µl **Wash Buffer**.
- Add 100 µl of **TMB Substrate Solution** to all wells including blank wells.
- Incubate the microwell strips for about 10-20 minutes at room temperature (18° to 25°C).
- Add 100 µl **Stop Solution** to all wells including blank wells.
- Blank microwell reader and measure colour intensity at 450 nm.

Note: Samples have been diluted 1:2, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 2).