

PRODUCT INFORMATION AND MANUAL

mouse TNF-RII Instant ELISA

BMS6003INST

Enzyme-linked immunosorbent assay for quantitative
detection of mouse TNF-RII.

For research use only.

Not for diagnostic or therapeutic procedures.

128 Tests

**mouse TNF-RII
BMS6003INST**



Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe

TABLE OF CONTENTS

1	Intended Use	3
2	Summary	3
3	Principles of the Test	4
4	Reagents Provided	5
5	Storage Instructions	5
6	Specimen Collection	6
7	Materials Required But Not Provided	7
8	Precautions for Use	8
9	Preparation of Reagents and Samples	10
10	Test Protocol	11
11	Calculation of Results	15
12	Limitations	18
13	Performance Characteristics	19
14	Bibliography	21
15	Ordering Information	24
16	Reagent Preparation Summary	25
17	Test Protocol Summary	26

1 Intended Use

The mouse TNF-RII Instant ELISA is an enzyme-linked immunosorbent assay for the quantitative detection of mouse soluble Type II TNF Receptor in cell culture supernatants, murine serum, plasma or other body fluids.. The mouse TNF-RII Instant ELISA is **for research use only. Not for diagnostic or therapeutic procedures.**

2 Summary

Tumor Necrosis Factor (TNF) was originally discovered in sera of animals and was found to cause hemorrhagic necrosis of some transplantable mouse and human tumors and to exhibit primarily cytotoxic activities against tumor but not normal cells in vitro (8, 14). The TNF family consists of two proteins designated TNF- α , also called cachectin (6), and TNF- β , also called lymphotoxin (15), which are pleiotropic cytokines that can mediate a wide variety of biological effects (3).

Both TNF- α and TNF- β have been shown to interact with a cell through specific high affinity receptors (4,8,9,12,16,17,18,21). TNF-receptors have been demonstrated on a wide variety of cells including fibroblasts (11), endothelial cells, adipocytes, liver membranes (5), granulocytes and several tumor cell lines (2,7,9,20). Normal and malignant myeloid cells as well as mitogen-stimulated lymphocytes express similar numbers of TNF receptors, whereas resting lymphoid cells have fewer, red blood cells and platelets have no detectable TNF receptors (13).

In most cases no correlation is observed between receptor number and sensitivity to TNF. In a variety of cell lines two different types of TNF receptors have been identified (7,10,19), TNF-RI and TNF-RII.

The present assay provides a simple, rapid and highly sensitive method for the determination of soluble mouse TNF-RII levels in body fluids or cell culture supernatants.

3 Principles of the Test

An anti-mouse TNF-RII monoclonal coating antibody is adsorbed onto microwells. Mouse TNF-RII present in the sample or standard binds to antibodies adsorbed to the microwells; a biotin-conjugated monoclonal anti-mouse TNF-RII antibody binds to mouse TNF-RII captured by the first antibody. Streptavidin-HRP binds to the biotin conjugated anti-mouse TNF-RII.

Following incubation unbound biotin conjugated anti mouse TNF-RII and Streptavidin-HRP is removed during a wash step, and substrate solution reactive with HRP is added to the wells.

A coloured product is formed in proportion to the amount of soluble mouse TNF-RII present in the sample. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from seven mouse TNF-RII standard dilutions and mouse TNF-RII sample concentration determined.

Figure 1

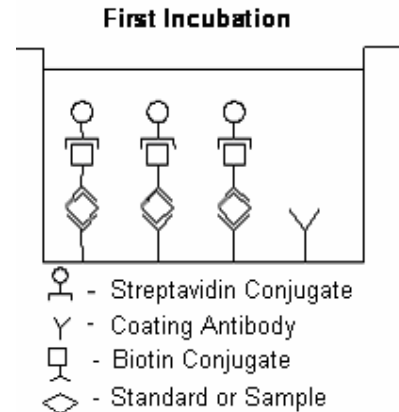


Figure 2


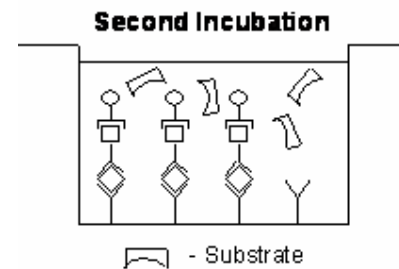


Figure 3

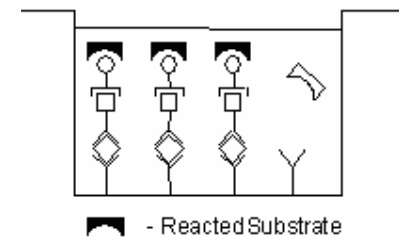

4 Reagents Provided

- 1 aluminium pouch with a **Microwell Plate coated with Monoclonal Antibody** (rat) to murine TNF-RII, **Biotin-Conjugate** (anti-mTNF-RII monoclonal antibody), **Streptavidin-HRP** and Sample Diluent, lyophilized
- 2 aluminium pouches with a **mouse TNF-RII Standard curve (coloured)**
- 1 bottle (25 ml) **Wash Buffer Concentrate** 20x (phosphate-buffered saline with 1% Tween 20)
- 1 vial (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 1 vial (50 ml) **Sample Diluent** (Use when an external predilution of the samples is needed)
- 1 vial (15 ml) **Stop Solution** (1M Phosphoric acid)
- 2 adhesive **Plate Covers**

5 Storage Instructions

Store ELISA plate and Standard curves whole kit at -20°C. The plate and the standard curves can also be removed, stored at -20°C, remaining kit reagents can be stored between 2° and 8°C. Expiry of the kit and reagents is stated on labels.

The expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, the reagent is not contaminated by the first handling.

6 Specimen Collection

Cell culture supernatants, murine serum, plasma or other biological samples are suitable for use in the assay. Remove the serum or plasma from the clot or red cells as soon as possible after clotting and separation.

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples must be stored frozen at -20°C to avoid loss of bioactive mouse TNF-RII. If samples are to be run within 24 hours, they may be stored at 2° to 8°C (for sample stability refer to 13).

Avoid repeated freeze-thaw cycles. Prior to assay, frozen serum or plasma should be brought to room temperature slowly and mixed gently.

7 Materials Required But Not Provided

- 5 ml and 10 ml graduated pipettes
- 5 μ l to 1000 μ l adjustable single channel micropipettes with disposable tips
- 50 μ l to 300 μ l adjustable multichannel micropipette with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform linear regression analysis

8 Precautions for Use

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statements(s) for specific advice.
- Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipette by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing substrate reagent.

- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5°C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

9 Preparation of Reagents and Samples

9.1 Wash Buffer

If crystals have formed in the Wash Buffer Concentrate, warm it gently until they have completely dissolved.

Pour entire contents (25 ml) of the Wash Buffer Concentrate into a clean 500 ml graduated cylinder. Bring to final volume to 500 ml with glass-distilled or deionized water. Mix gently to avoid foaming. The pH of the final solution should adjust to 7.4.

Transfer to a clean wash bottle and store at 2° to 25°C. Please note that Wash Buffer is stable for 30 days.

10 Test Protocol

- **Use plate immediately after removal from -20°C!**
 - **Do not wait until pellets have completely dissolved before applying samples - the binding reaction in the standard strips starts immediately after addition of water!**
 - **Do not try to dissolve pellets by pipetting up and down in the wells - some parts of the pellet could stick to the tip creating high variation of results**
 - **Perform the washing step with at least 400 µl of washing buffer as stated in the manual or fill the wells completely - otherwise any pellet residues sticking to the rim of the well will not be removed and create high variation of results**
 - **Allow the washing buffer to sit in the wells for a few seconds before aspiration**
 - **Remove covers of the standard strips carefully in order that all the lyophilised pellets remain in the wells**
- a. Prepare your samples before starting with the test procedure. Dilute serum or plasma samples 1:50 with Sample Diluent according to the following dilution scheme:
5 µl Sample + 245 µl Sample Diluent
 - b. Determine the number of microwell Strips required to test the desired number of samples plus microwell Strips for blanks and standards (coloured). Each sample, standard, blank, and optional control sample should be assayed in duplicate. Remove extra microwell Strips from holder and store in foil bag with the desiccant provided at -20°C sealed tightly. Place microwell strips containing the standard curve in position A1/A2 to H1/H2 (see Table 1).
 - c. Add 100 µl of distilled water to the sample wells.
 - d. Add distilled water to all standard and blank wells as indicated on the label of the standard strips (A1, A2 to H1, H2).

Table 1

Table depicting an example of the arrangement of blanks, standards and samples in the microwell strips:

	1	2	3	4
A	Standard 1 (500 pg/ml)	Standard 1 (500 pg/ml)	Sample 1	Sample 1
B	Standard 2 (250 pg/ml)	Standard 2 (250 pg/ml)	Sample 2	Sample 2
C	Standard 3 (125 pg/ml)	Standard 3 (125 pg/ml)	Sample 3	Sample 3
D	Standard 4 (62.5 pg/ml)	Standard 4 (62.5 pg/ml)	Sample 4	Sample 4
E	Standard 5 (31.25 pg/ml)	Standard 5 (31.25 pg/ml)	Sample 5	Sample 5
F	Standard 6 (15.6 pg/ml)	Standard 6 (15.6 pg/ml)	Sample 6	Sample 6
G	Standard 7 (7.8 pg/ml)	Standard 7 (7.8 pg/ml)	Sample 7	Sample 7
H	Blank	Blank	Sample 8	Sample 8

- e. Add 50 µl of **1:50 prediluted Sample**, in duplicate, to the designated wells and mix the contents.
- f. Cover with a **Plate Cover** and incubate at room temperature (18°C to 25°C) for 3 hours, if available on a microplate shaker at 200 rpm.
- g. Remove **Plate Cover** and empty wells. Wash the microwell strips 3 times with approximately 400 µl Wash Buffer per well with thorough aspiration of microwell contents between washes. Take care not to scratch the surface of the microwells.

After the last wash, tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing or place upside down on a wet absorbent paper for no longer than 15 minutes. Do not allow wells to dry.

- h. Pipette 100 µl of **TMB Substrate Solution** to all wells, including the blank wells.
- i. Incubate the microwell strips at room temperature (18° to 25 °C) for about 10 min. Avoid direct exposure to intense light.

The colour development on the plate should be monitored and the substrate reaction stopped (see point j. of this protocol) before positive wells are no longer properly recordable. Determination of the ideal time period for colour development has to be done individually for each assay.

It is recommended to add the stop solution when the highest standard has developed a dark blue colour. The colour development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.6 – 0.65.

- j. Stop the enzyme reaction by quickly pipetting 100 µl of **Stop Solution** into each well, including the blank wells. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2 - 8°C in the dark.
- k. Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the

reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the mouse TNF-RII standards.

Note: In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.

11 Calculation of Results

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20 per cent of the mean.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the mouse TNF-RII concentration on the abscissa. Draw a best fit curve through the points of the graph.
- To determine the concentration of circulating mouse TNF-RII for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding mouse TNF-RII concentration.
- ***Samples have been diluted 1:100, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 100).**
- It is suggested that each testing facility establishes a control sample of known mouse TNF-RII concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.
- A representative standard curve is shown in Figure 4. This curve cannot be used to derive test results. Every laboratory must prepare a standard curve for each group of microwell strips assayed.

* N.B: There is a common dilution factor for samples due to the conjugate which must then be included in the calculation. The samples contribute 100 µl to the final volume per well. These 100 µl are composed of 50 µl of sample diluent plus 50 µl of the 1:50 prediluted sample. This is a 1:100 dilution.

The remaining 50 µl to give 150 µl are due to the addition of 50 µl conjugate to all wells.

50 µl sample diluent and 50 µl conjugate results in 100 µl reconstitution volume, addition of 50 µl of the 1:50 prediluted sample (50 µl + 50 µl of 1:50 prediluted sample = 1:100 dilution)

Figure 4

Representative standard curve for mouse TNF-RII Instant ELISA. Mouse TNF-RII was diluted in serial 2-fold steps in Sample Diluent, each symbol represents the mean of 3 parallel titrations. Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

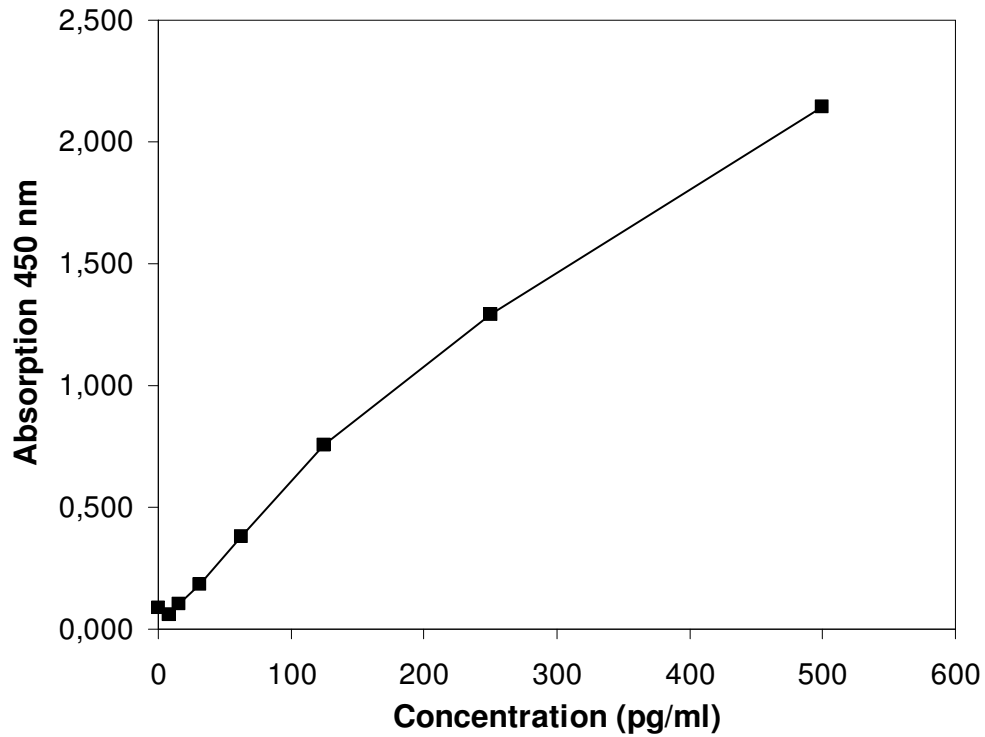


Table 2

Typical data using the mouse TNF-RII INSTANT ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	Mouse TNF-RII Concentration (pg/ml)	O.D. (450 nm)	O.D. Mean	C.V. (%)
1	500	2.102	2.153	3.3
	500	2.204		
2	250	1.328	1.301	3.0
	250	1.273		
3	125	0.730	0.766	6.6
	125	0.802		
4	62.5	0.396	0.388	3.1
	62.5	0.379		
5	31.25	0.179	0.191	8.9
	31.25	0.203		
6	15.6	0.110	0.113	3.1
	15.6	0.115		
7	7.8	0.069	0.068	2.1
	7.8	0.067		
Blank		0.008	0.009	
		0.009		

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus colour intensity. Values measured are still valid.

12 Limitations

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.
- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.

13 Performance Characteristics

13.1 Sensitivity

The limit of detection of mouse TNF-RII defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be 0.6 pg/ml (mean of 6 independent assays).

13.2 Reproducibility

13.2.1 Intra-assay

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 serum samples containing different concentrations of mouse TNF-RII. 2 standard curves were run on each plate. The calculated overall intra-assay coefficient of variation was <10%.

13.2.2 Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments by 3 technicians. Each assay was carried out with 6 replicates of 8 serum samples containing different concentrations of mouse TNF-RII. 2 standard curves were run on each plate. The calculated overall inter-assay coefficient of variation was <10%.

13.3 Spike Recovery

The spike recovery was evaluated by spiking 4 levels of mouse TNF-RII into normal murine serum. Recoveries were determined in 3 independent experiments with 4 replicates each. The unspiked serum was used as blank in these experiments. The overall mean recovery was 90.3%.

13.4 Dilution Parallelism

4 serum samples with different levels of mouse TNF-RII were analysed at serial 2 fold dilutions with 4 replicates each. The overall mean recovery was 113.3%.

13.5 Sample Stability

13.5.1 Freeze-Thaw Stability

Aliquots of serum samples (unspiked or spiked) were stored at -20 °C and thawed 5 times, and the mouse TNF-RII levels determined. There was no significant loss of mouse TNF-RII immunoreactivity detected by freezing and thawing.

13.5.2 Storage Stability

Aliquots of serum samples (spiked or unspiked) were stored at -20 °C, 2-8 °C, room temperature (RT) and at 37 °C, and the mouse TNF-RII level determined after 24 h. There was no significant loss of mouse TNF-RII immunoreactivity detected during storage at -20 °C, 4 °C and room temperature. Storage at 37 °C gave rise to about 50 % loss of mTNF-RII immunoreactivity.

13.6 Specificity

To define the specificity of this ELISA several proteins were tested for cross reactivity. There was no cross reactivity observed.

14 Bibliography

- 1) Adolf G., I. Apfler. (1991). A monoclonal antibody-based enzyme immunoassay for quantitation of human tumor necrosis factor binding protein I, a soluble fragment of the 60 kDa TNF receptor, in biological fluids. *J. Immunol. Meth.* 143, 127-136.
- 2) Aggarwal B., T. Eessalu, and P. Hass. (1985). Characterization of receptors for human tumour necrosis factor and their regulation by gamma-interferon. *Nature* 318, 665-667.
- 3) Aggarwal B. (1987). Tumour necrosis factor - TNF- α . and TNF- β : their structure pleiotropic biological effects. *Drugs Future* 12, 891-898.
- 4) Baglioni C., S. McCandless, J. Tavernier, and W. Fiers (1985). Binding of human tumor necrosis factor to high affinity receptors on HeLa and lymphoblastoid cells sensitive to growth inhibition. *J. Biol. Chem.* 260, 13395-13397.
- 5) Beutler B., J. Mahoney, N. Le Trang, P. Pekala, and A. Cerami. (1985). Purification of cachectin, a lipoprotein lipase-suppressing hormone secreted by endotoxin-induced RAW 264.7 cells. *J. Exp. Med.* 161, 984-995.
- 6) Beutler B., and A. Cerami. (1987). Cachectin: more than a tumor necrosis factor. *N. Engl. J. Med.* 316, 379-385.
- 7) Brockhaus M., H. Schoenfeld, E. Schlaeger, W. Hunziker, W. Lesslauer, and H. Loetscher. (1990). Identification of two types of tumor necrosis factor receptors on human cell lines by monoclonal antibodies. *Proc. Natl. Acad. Sci. USA* 87, 3127-3131.
- 8) Carswell E., L. Old, R. Kassel, S. Green, N. Fiore, and B. Williamson. (1975). An endotoxin-induced serum factor that causes necrosis of tumors. *Proc. Natl. Acad. Sci. USA* 72, 3666-3670.
- 9) Creasey A., R. Yamamoto, and C. Vitt. (1987). A high molecular weight component of the human tumor necrosis factor receptor is associated with cytotoxicity. *Proc. Natl. Acad. Sci. USA* 84, 3293-3297.

- 10) Hohmann H., R. Remy, M. Brockhaus, and A. van Loon. (1989). Two different cell types have different major receptors for human tumor necrosis factor (TNF- α). *J. Biol. Chem.* 264, 14927-14934.
- 11) Kull F., S. Jacobs, and P. Cuatrecasas. (1985). Cellular receptor for 125I-labeled tumor necrosis factor: specific binding, affinity labeling, and relationship to sensitivity. *Proc. Natl. Acad. Sci. USA* 82, 5756-5760.
- 12) Loetscher H., Y. Pan, H.-W. Lahm, R. Gentz, M. Brockhaus, H. Tabuchi, and W. Lesslauer. (1990). Molecular Cloning and Expression of the Human 55 kd Tumor Necrosis Factor Receptor. *Cell* 61, 351-359.
- 13) Munker R., J. DiPersio, and H. Koeffler. (1987). Tumor Necrosis Factor: receptors on hematopoietic cells. *Blood* 70, 1730-1734.
- 14) Old L. (1985). Tumor necrosis factor (TNF). *Science* 230, 630-632.
- 15) Paul N., and N. Ruddle. (1988). Lymphotoxin. *Ann. Rev. Immunol.* 6, 407 - 438.
- 16) Rubin B., S. Anderson, S. Sullivan, B. Williamson, E. Carswell, and L. Old. (1985). High affinity binding of 125I-labeled human tumor necrosis factor (LuKII) to specific cell surface receptors. *J. Exp. Med.* 162, 1099-1104.
- 17) Scheurich P., U. Ücer, M. Krönke, and K. Pfizenmaier. (1986). Quantification and characterization of high-affinity membrane receptors for tumor necrosis factor on human leukemic cell lines. *Int. J. Cancer* 38, 127-133.
- 18) Scheurich P., B. Thoma, U. Ücer, and K. Pfizenmaier. (1987). Immunoregulatory activity of recombinant human tumor necrosis factor (TNF)- α .: Induction of TNF receptors on human T cells and TNF- α . mediated enhancement of T cell responses. *J. Immunol.* 138, 1786-1790.
- 19) Sprang S. (1990). The divergent receptors for TNF. *TIBS* 15, 366-368.

- 20) Thoma B., M. Grell, K. Pfizenmaier, and P. Scheurich. (1990). Identification of a 60-kD tumor necrosis factor (TNF) receptor as the major signal transducing component in TNF responses. *J. Exp. Med.* 172, 1019-1023.
- 21) Tsujimoto M., Y. Yip, and J. Vilcek. (1985). Tumor necrosis factor: Specific binding and internalization in sensitive and resistant cells. *Proc. Natl. Acad. Sci. USA* 82, 7626-7630

For literature update refer to **www.bendermedsystems.com/44.html**

15 Ordering Information

For orders please contact:

Europe-Headquarters

Bender MedSystems GmbH
Campus Vienna Biocenter 2
A-1030 Vienna, Austria, Europe
phone: +43 1 796 40 40 ext. 114
fax: +43 1 796 40 40 ext. 400
e-mail: order@bendermedsystems.com

USA

Bender MedSystems, Inc.
849 Hinckley Road
Burlingame, CA 94010, USA
toll-free phone: +1 (866) 952 2112
toll-free fax: +1 (877) 952 2112
phone: +1 (650) 952 2112
fax: +1 (650) 952 2252

For technical information please contact:

e-mail: techserv@bendermedsystems.com
www.bendermedsystems.com

Cat.No. BMS6003INST mouse TNF-RII INSTANT ELISA

16 Reagent Preparation Summary

16.1 Wash Buffer

Add **Wash Buffer Concentrate** 20 x (25 ml) to 475 ml distilled water

17 Test Protocol Summary

- Predilute sample with Sample Diluent 1:50.
- Place standard strips in position A1/A2 to H1/H2.
- Add 100 µl **distilled water** to sample wells.
- Add **distilled water**, in duplicate, to all standard and blank wells as indicated on the label of the standard strips.
- Add 50 µl **Sample** to designated wells.
- Cover microwell strips and incubate 3 hours at room temperature (18° to 25°C) on a microplate shaker at 200 rpm.
- Empty and wash microwell strips 3 times with 400 µl **Wash Buffer**.
- Add 100 µl of **TMB Substrate Solution** to all wells including blank wells.
- Incubate the microwell strips for about 10 minutes at room temperature (18° to 25°C).
- Add 100 µl **Stop Solution** to all wells including blank wells.
- Blank microwell reader and measure colour intensity at 450 nm.

Note: Samples have been diluted 1:100, thus the concentration read from the standard curve must be multiplied by the dilution factor (x 100).